




Chemistry



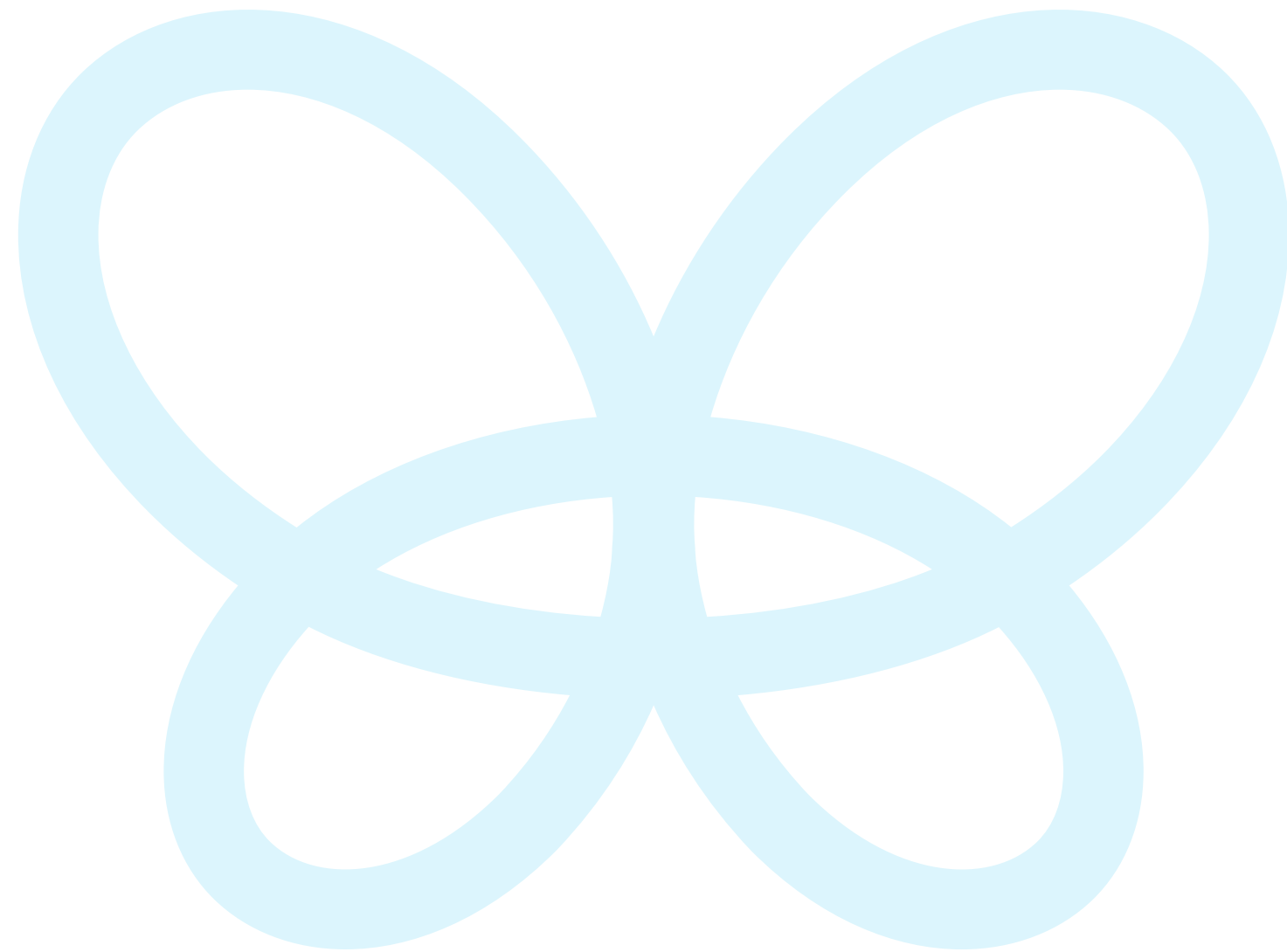

Flexible and ductless filtration solutions

to ensure safety in laboratories
and enclosed public spaces

Who are we?

Since 1968, **Erlab** has been the world-leading specialist and inventor of **standalone, ductless and zero-emission filtering fume hoods**, enabling chemicals to be handled safely in laboratories.

Since day one, we have been entirely focused on researching, designing, developing and manufacturing long-lasting protection solutions. Our primary objective is to provide our users with the most effective solutions to protect them against chemical inhalation hazards in laboratories. Our worldwide presence, production facilities and strong research and development activity enable us to offer the most advanced filtration technology to laboratories in the chemical, pharmaceutical, cosmetics and food industries, as well as to hospitals and educational institutions.



Left: Stéphane Hauville: Chairman and CEO
Right: Antoine Hauville: Managing Director
Centre: Marie Hauville: Co-founder of Erlab

AFNOR* NF X 15-211:2009 standard

Compliance with standards is of utmost importance to us. Based on scientific criteria, the **AFNOR NF X 15-211:2009** reflects the level of performance of our products, and ensures you are protected every day at your workstation.

Our expertise in recirculating filtered air places the environment at the heart of your laboratory: all our solutions have been designed to limit environmental impact and contribute towards one of today's major global challenges: saving energy.

*French Standardization Association [Association Française de Normalisation]



François P. Hauville: 1932 - 2011: Founder of Erlab

International deployment from day one.

Filtering toxic gases and hazardous particles is a global issue.

We have installed over **150 000** ductless filtering fume hoods in laboratories around the world.



Europe: Erlab S.A.S. (France)



North America: Erlab, Inc. (USA)



Asia: Erlab Ltd (China)



Our three Erlab production sites: France, United States and China




ERLAB ABOVE Filtration Technology

The result of over 50 years of research and innovation, **ERLAB ABOVE FILTRATION TECHNOLOGY** is a mark of quality and a safety guarantee in the field of filtration technology which is **used to protect laboratory personnel**.

We use our technological developments and expertise to drive solutions, transforming the impossible into the possible.

ERLAB ABOVE is the invisible difference. It brings together technology that filters, detects and communicates, making the laboratory air you breathe cleaner and safer, both inside and out.

Accept nothing less than **ERLAB ABOVE** filtration technology. A proven mark of quality since 1972.

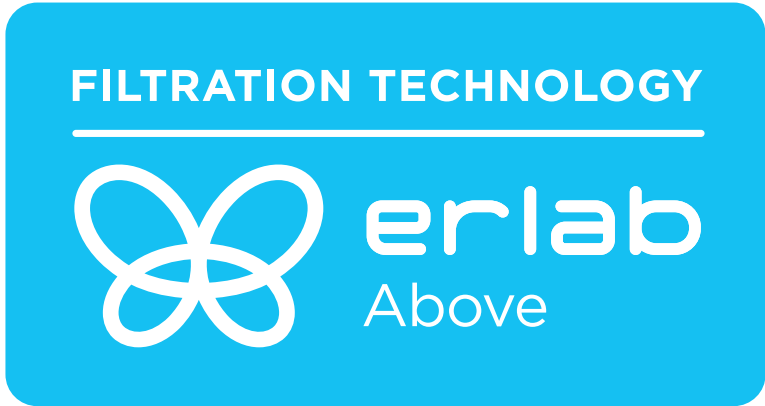

CONTENTS

THE BENEFITS OF ERLAB PRODUCTS..... 08

- Safety
- Protect the environment
- Reduce installation costs
- Deliver energy savings
- Optimise the flexibility and design of laboratories
- Proven technology

THE QUALITY OF ERLAB PRODUCT DESIGN..... 10

- Research and development
- Design
- Production
- Guarantees

THE ESP..... 12

STANDARDS..... 14

- Filtration efficiency
- Containment efficiency of the enclosure

FILTRATION AND FILTRATION TECHNOLOGY..... 18

DETECTION..... 22

FLEX TECHNOLOGY 24

SMART TECHNOLOGY 26

ERLAB PRODUCTS AND THEIR APPLICATIONS..... 28

Chemical handlings

- Captair filtering fume hoods 30
- Captair weighing stations 40

HANDLINGS AND SAMPLES PROTECTION

- Biological hoods 48
- Filtering clean air workstations..... 54

Chemical storage

- Filtering storage cabinets 62
- Filtration boxes..... 72

Lab air pollution

- Halo air purifiers..... 78
- Halo Sense air quality sensor 86

Complete protection of investigations

- The Pyramid glove bag 92

SERVICES 96

- eGuard
- eValiQuest
- Erlab Maintenance
- Erlab filters

THE BENEFITS OF ERLAB PRODUCTS

Erlab solutions are designed to ensure people's safety when handling chemicals in the laboratory. Based on the principle of filtration, they offer increased protection against the risk of chemical inhalation caused by the emission of harmful molecules and particles at the workstation.

Safety

Erlab solutions are designed to ensure people's safety when handling chemicals in the laboratory. Based on the principle of filtration, they offer Erlab increased protection against the risk of chemical inhalation caused by the emission of harmful molecules and particles at the workstation.



Protect the environment

Erlab solutions do not require a connected ventilation system, which stops pollutants from being released straight into the atmosphere and therefore helps protect the environment. Furthermore, they do not create pollution like traditional extraction hoods which need to produce energy to keep the ventilation systems functioning correctly.

Reduce installation costs

Erlab solutions are quick and easy to install. Unlike connected systems, they do not require a ventilation system linked to an air supply and extraction system to be installed.

You just need one plug to get up and running, and you can get set up at any suitable moment, without complex planning.

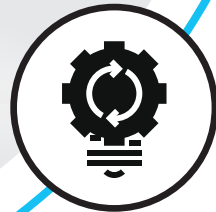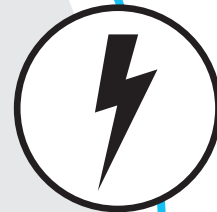

Deliver energy savings

Creating ventilation balance is essential for connected systems to function correctly; however, this requires significant energy consumption. Erlab solutions, on the other hand, incur no energy costs for the extraction and supply of conditioned air, nor do they have high operational costs, even taking into account the cost to change the filter.

Optimise the flexibility and design of laboratories


Erlab solutions can be easily moved within a laboratory to suit your protection needs, and without affecting the ventilation balance of the room.

Proven technology

With over 150,000 protection enclosures installed around the world, Erlab offers a proven product ecosystem.

THE QUALITY OF ERLAB PRODUCT DESIGN

Our products are the result of 50 years of research and development. Your safety is our primary concern.

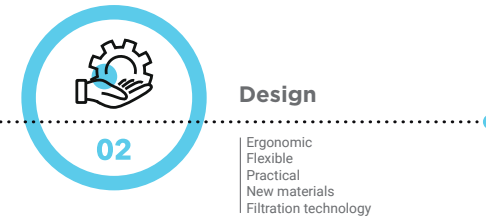

Research and development

Our goal is to be, and to continue to be, the world's leading partner in personal safety. At Erlab, we are therefore continually investing in our state-of-the-art laboratories and research teams. All of our Erlab engineers and chemists are experts in molecular filtration. In the Erlab laboratory, they use sophisticated analytical equipment to conduct the hundreds of tests required to meet the standards, as well as to determine and optimise the retention capacity of the filters. We continuously monitor the safety of our Erlab products, constantly improve the quality and technologies we used, as well as develop new products that offer enhanced protection for laboratory personnel.

Equipment for filtration and containment tests worth more than 1 million.



Design

At Erlab, we place great emphasis on the design of our products. For this reason, the dimensions, visibility, front openings, mounting and working position are designed to best meet the needs of those who handle samples.

Enclosures

- Enclosure width: from 80 cm to 180 cm
- Enclosure depth: from 621 cm to 960 cm
- Easy to install large instruments
- Easy integration with laboratory furniture

Visibility

- High optical quality synthetic glass panels for an exceptional view of handlings
- Lighting suitable for precise work conditions
- Tilted front panel for excellent visibility

Front panel openings

- Large range of movement in the enclosure
- Central protection screen avoids any risk of chemical projection

Installation

- Ready to install, quick to assemble
- Very few tools required

Working position

- Fatigue-free handling while sitting or standing
- Work surface with rounded edges to provide rest for forearms
- Tilted front panel provides visual comfort and a comfortable

Home machine interface

- Simple, intuitive and connectable
- Monitor safety in real time



Production

Production has always been at the heart of Erlab's operations. This is an essential area, highlighted by the work of the design team, who produce and optimise our solutions using the most advanced production technologies, in order to meet market expectations.



Guarantees and commitments

Erlab only provides protective enclosures to customers once they have verified that the enclosure is suitable for the chemicals they wish to use. The customer's specific handlings are first analysed in Erlab's internal testing laboratory to determine if they can be carried out safely with a filtering fume hood from the Erlab range, before Erlab equips their laboratory.

THE ESP®

A package of three services is included with your purchase, designed to ensure your safety throughout the service life of your device.


WHAT IS THE ERLAB SAFETY PROGRAM?

Erlab's long-term commitment to ensuring your safety.

Erlab's R&D laboratory analyses the interactions between molecules and particles in order to approve the filtration technology for your handlings. Based on this scientific analysis, our laboratory will recommend the type of device, filtration column configuration and enclosure size that will ensure your protection.

Contact your **ESP** specialist today to set up your Erlab protection solution.

Website: www.erlab.com

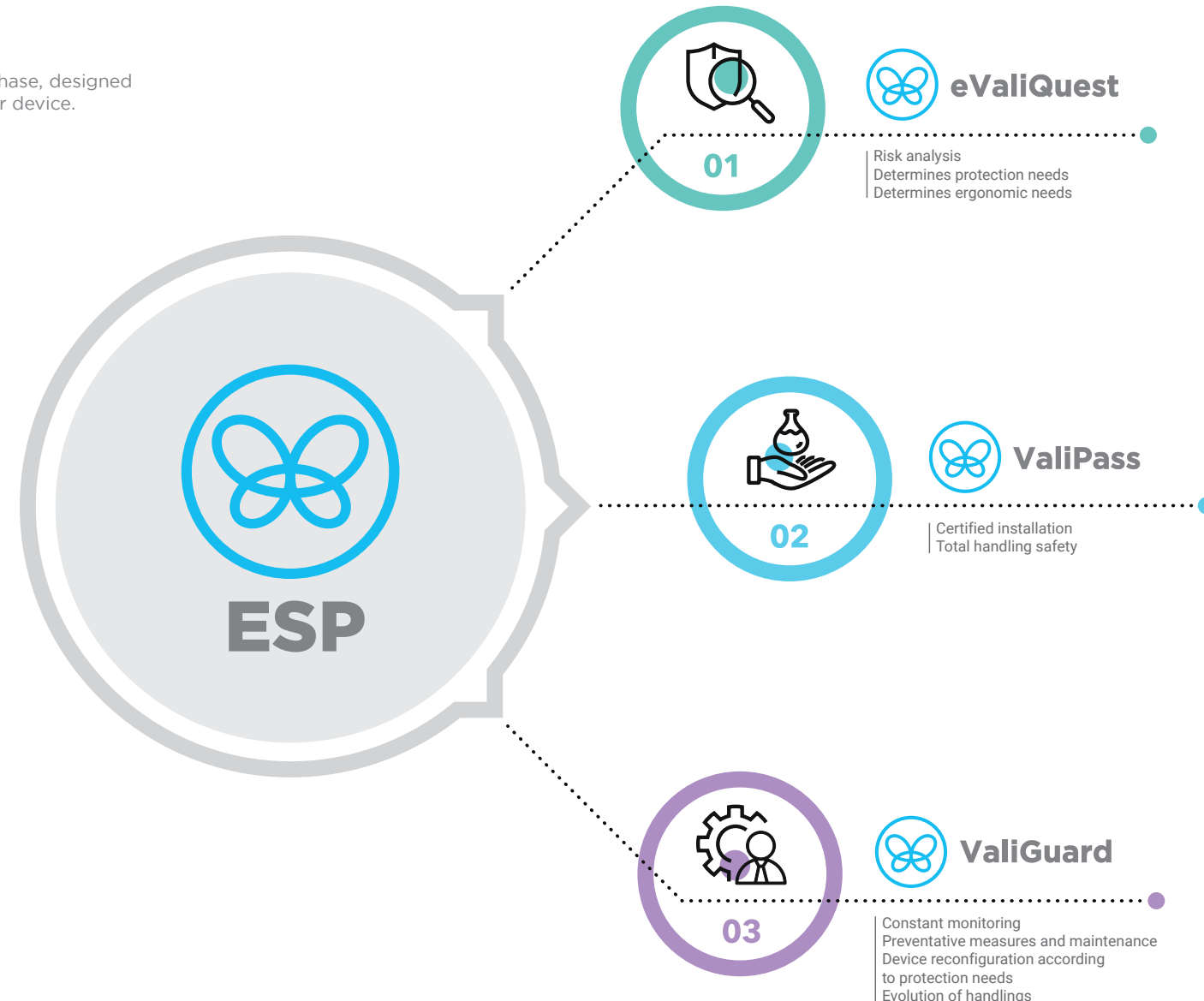

THE ESP®




eValiQuest®

An ESP agent will help you fill out the initial questionnaire where you will specify the handlings that you wish to carry out. Within 48 hours, our laboratory specialists will suggest which type of device and filtration technology will be suitable for your use. We are committed to ensuring your protection by certifying the feasibility of your handlings.



ValiPass®

When your fume hood is installed, a usage certificate will indicate the specific chemicals used, the type of filter and the estimated service life for which your device has been validated. This certificate acts as a permanent reminder for the user or the person responsible for data safety for the use of the device.

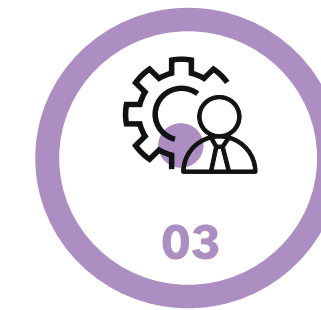

ValiGuard®

An ESP® agent will contact you periodically to ensure that your handlings have not changed and that the filter is still effective. They will guide you through each step for checking the performance of the filter as well as for replacing it. If a change to handlings is noted, the ESP® agent will ask you to complete a new questionnaire (see step 1). After the assessment, a new usage certificate specifying the authorised chemical products will be sent to you, in order to carry out handlings in optimal safety conditions.

STANDARDS

The AFNOR NF X 15-211 standard

This is currently the most advanced and rigorous industry standard which assesses the safety of filtering fume hoods. The French Union of Mechanical Standardisation, comprising a panel of experts (from the French National Research and Safety Institute, state bodies and professional unions), was appointed by **AFNOR** to establish the AFNOR NF X 15-211:2009 standard. This standard applies to filtering fume hoods (also known as recirculating fume cupboards or **Enclosures for Toxics using Recirculating Air Filtration**) which are designed for research work, analysis, teaching, etc. in any laboratory in which chemical agents subject to Threshold Limit Values (TLV) or Occupational Exposure Limits (OEL) are handled.

This standard sets out performance requirements related to:



- Filtration efficiency
- Containment efficiency
- Air face velocity
- Documentation: **CHEMICAL LISTING***

** Guide to filtered or retained chemicals*

The classes set out by the standard are:


Class 1: Filtering fume hood with safety reserve, one main filtration level, one safety filtration level and continuous monitoring of filtration efficiency.

Class 2: Filtering fume hood without safety reserve and with one main filtration level.

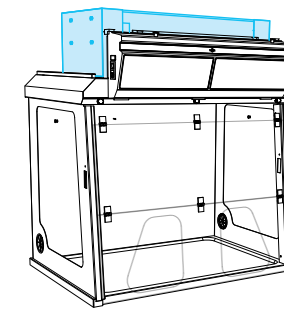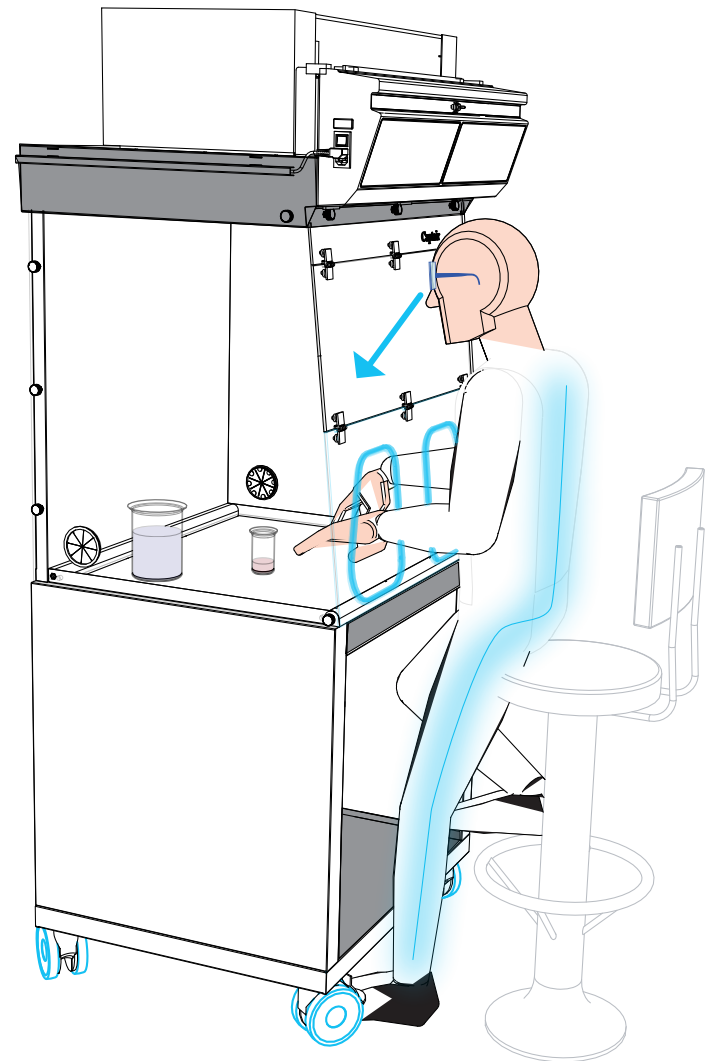


Filtration efficiency

This is defined as the filter's ability to trap harmful molecules handled in the enclosure and determines the quality of the recirculated air downstream of the filter.

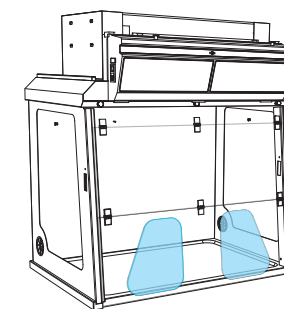


Air face velocity

This refers to the ability of the fume hood to create a dynamic air barrier between the operator and the handling.

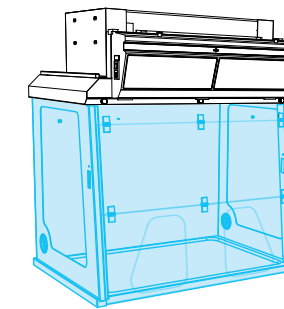

Containment efficiency

This is defined as the capacity of the fume hood to keep vapours or particles within the enclosure, preventing their release into the laboratory environment.

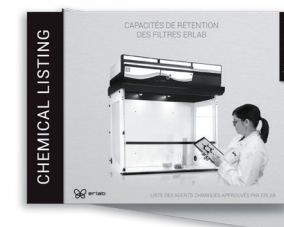

Documentation: Chemical Listing*

Erlab has created a **CHEMICAL LISTING** a guide to over 700 chemicals. This guide is provided with every enclosure, in compliance with the standard.


** Guide to filtered or retained chemicals*



Filtration efficiency


Erlab’s filtration technology conforms with the **NF X 15-211:2009 standard**, the most rigorous industry standard for molecular filtration, developed by a committee of independent scientists and specialised manufacturers.

This standard establishes the performance criteria which impose a maximum release of 1% of the TLV (Threshold Limit Value) of the handled products.



Solvents test

Performance criterion met: <1% of the TLV



Acids test

Performance criterion met: <1% of the TLV






Class 1		Class 2
Normal operating phase	Normal operating phase, during which the concentration downstream of the filters is less than 1% of the TLV	
Detection phase	Detection phase, during which the concentration downstream of the filters must be less than 1% of the TLV and during which the automatic filter failure detector must alert the user	Detection phase, during which the concentration downstream of the filters must be less than 50% of the TLV
Safety phase	Safety phase, during which the concentration downstream of the filters must be less than 50% of the TLV and whose duration must not be less than 1/12 of the duration of the normal operating phase	

Containment efficiency of the enclosure

The level of containment is defined by the capacity of the fume cupboard to keep pollutants within the enclosure, preventing their release into the laboratory environment.

Containment efficiency is proven by tests carried out according to the protocols outlined in the **EN 14175-3, ASHRAE 110:2016 and AFNOR NF X 15-211:2009 standards**.

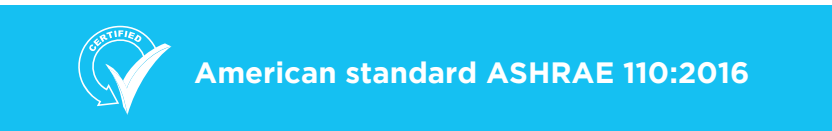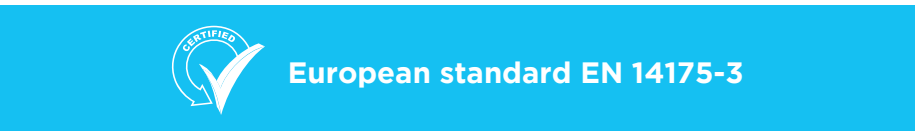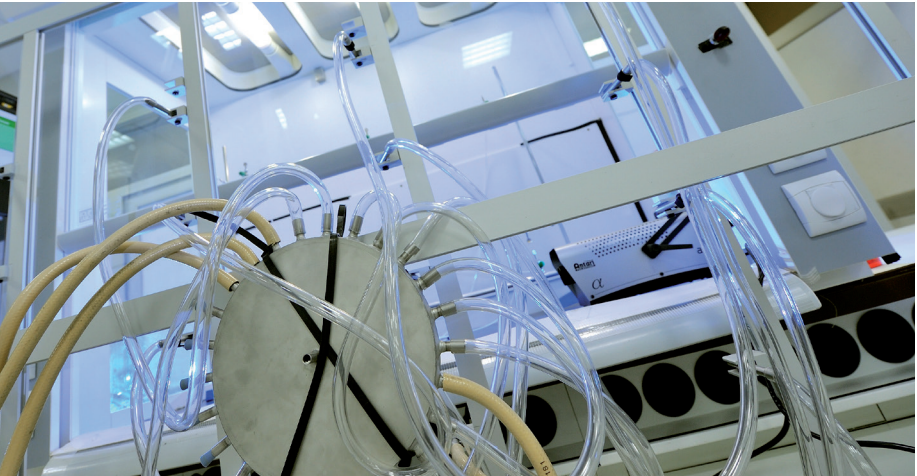

FILTRATION

Expertise and technologies developed by a leading R&D team.

We provide filtration technology that offers a high level of protection for laboratory personnel from inhaling chemical substances.



Our unique solutions allow pollution to be captured at source and in the filters before it is released into the clean air of the work environment. These filter cartridges are available as part of a unique activated carbon range, which is designed to protect personnel from inhalation.


This is possible thanks to the filtration technology that our Research and Development department has been continuously improving for over 50 years.

A very strict product specification based on compliance with standardised international protocols allows us to select raw materials and develop technologies with tailored porosity which have the capacity, under normal usage conditions, to adsorb a very broad spectrum of molecules without risk of desorption, just like in military-style gas masks. Our filters are subject to rigorous tests in compliance with the **AFNOR NF X 15-211:2009 standard**, the reference in the field of filtering fume hoods.


?

What is adsorption?


Adsorption is a group of physico-chemical surface reactions through which free molecules are condensed on the surface of a solid. They can take place in a liquid or gaseous medium. The phenomenon of adsorption has long been of interest for the capture of gaseous pollutants, particularly in gas masks or filtering fume cupboards. Be careful, however, not to confuse adsorption with absorption.

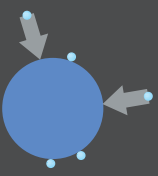

ADSORPTION

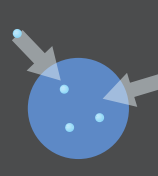

ABSORPTION



FILTRATION TECHNOLOGY

Carbon filter anti-detachment system US patent number 7,563,301.

Air pollutants in your laboratory

Chemical agents, in the form of gases and/or particles, are an inhalation health risk to those who work in laboratories. The health authorities have established concentration thresholds that must not be exceeded under any circumstances and that are defined by occupational exposure limit values (TLV), expressed in parts per million (PPM). The harmful omnipresence of these pollutants due to everyday handling requires all laboratories to adopt prevention and protection measures which comply with the regulations in force.

With over 50 years of expertise in filtration technology, **Erlab** has developed Flex technology, which combines molecular and particulate filtration technologies to provide a global protection solution for the most common handlings encountered in all laboratory disciplines, regardless of their environment and sector of activity.

Molecular filtration technology: super-activated carbon

For over a century, activated carbon has been used for its exceptional adsorption properties. The different varieties are used today in multiple applications, including water treatment, VOC, solvent recovery and chemical catalysis.

Each of these uses requires an activated carbon with unique and adapted physical and chemical properties.

Our experience is based on more than 50 years of testing and is expressed in our Chemical Listing guide, which attests to our excellent mastery of filtration.

We also integrate an environmental dimension when we develop our filtration technology; in particular, refusing to use environmentally toxic substances, which we have avoided for many years.

Our filters are subject to rigorous tests in compliance with the **AFNOR NF X 15-211:2009**. Their performance levels, measured by the results obtained, are a guarantee of safety for those who use our solutions.

In terms of safety, each of our filters is delivered with a quality certificate that tracks its entire manufacturing cycle.

Which molecular filter for which handling?

Different filter types	
AS	For organic vapours
BE/BE+	Multipurpose for acid vapours + organic vapours
F	For formaldehyde vapours
K	For ammonia vapours


Particulate filtration technology: HEPA 14

This filtration technology traps particles larger than 0.1 μm in diameter with **an efficiency of 99.995% according to the MPPS (Most Penetrating Particle Size) method of the EN 1822-1 standard.**

ERLAB FILTERS


Erlab filters contain the world's most powerful carbon.

Only Erlab manufactures high efficiency filters, the result of 50 years of research and development.

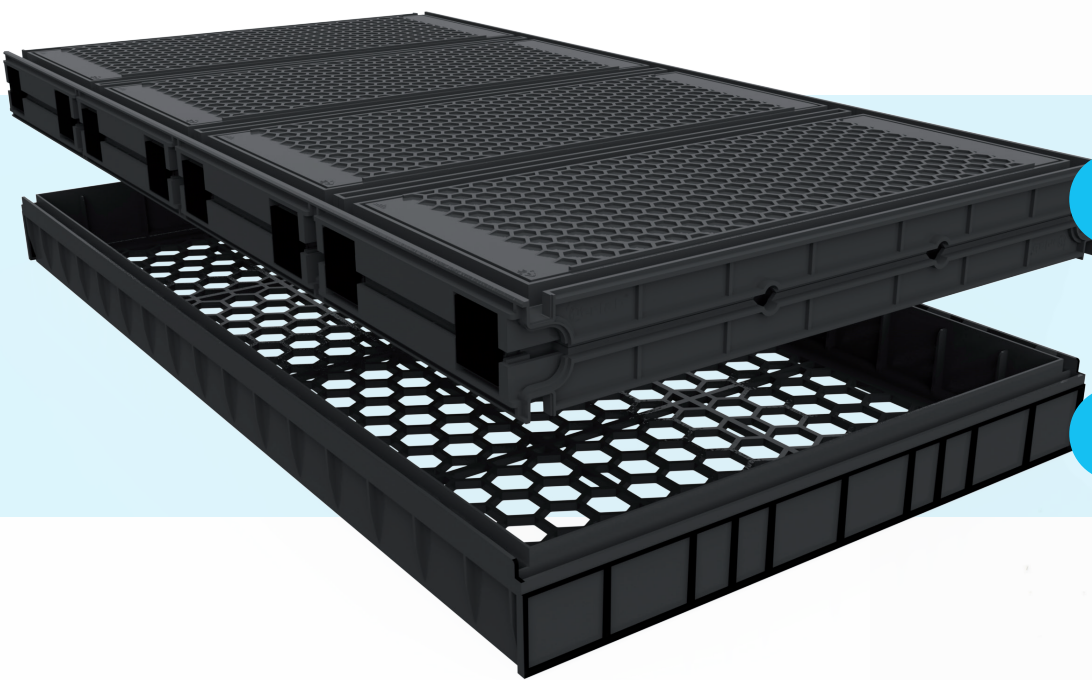

1 Activated carbon filtration zone. Anti-pollutant case with carbon bed height control.

2 Prefiltration zone.

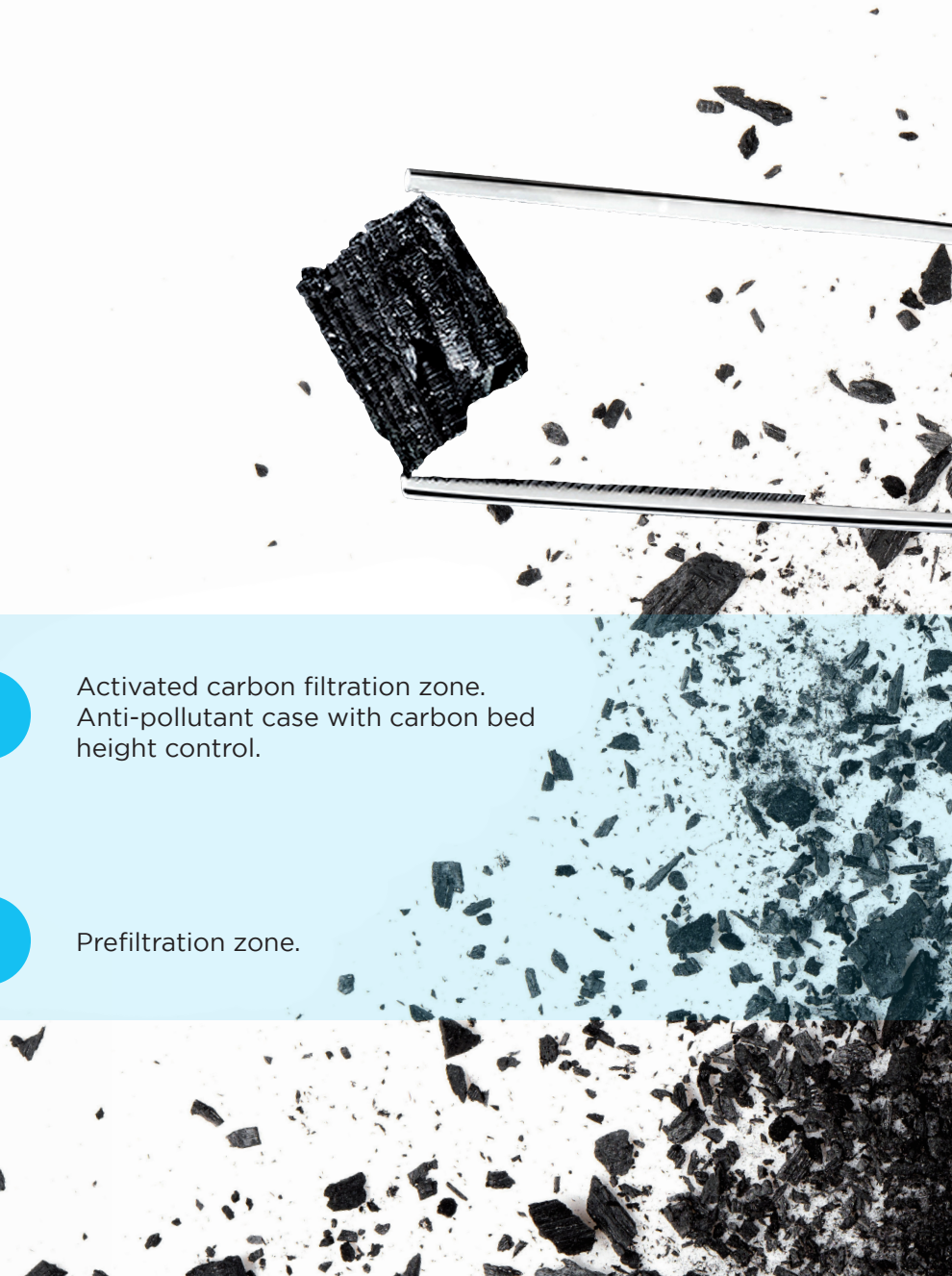

REVOLVING® FILTER SYSTEM

Optimising the main filters' service life

The advantages of the Revolving system

Significant optimisation of the main filter's service life
Substantial savings on replacement costs.



The patented concept Revolving filter system

“Class 1” filtration column has two filtration levels:

One main module + one safety module.

The Revolving system allows you to replace one filter module at a time and optimises the main filter's service life. No pollutants are released into the laboratory, even if there is a fault with the main filter.

When the main filter approaches saturation, the few non-retained molecules are directed straight towards the safety filter. The safety filter replaces the main filter when the latter reaches maximum capacity. A new filter is then installed in place of the safety filter.

DETECTION


Sensors

Erlab's exclusive detection system consists of three sensors which detect very broad spectrum of molecules. This allows users to monitor filtration efficiency and performance.

The three sensors: **VOC** (Molecode S), **Formaldehyde** (Molecode F) and **Acids** (Molecode A) are integrated following the automatic online chemical hazard analysis via **our eValiQuest service**.

The detection system also has other sensors which check your enclose is functioning correctly, in real time.

You can find this information on the device's embedded service (p.26-27) and on the eGuard® management interface (p.94).

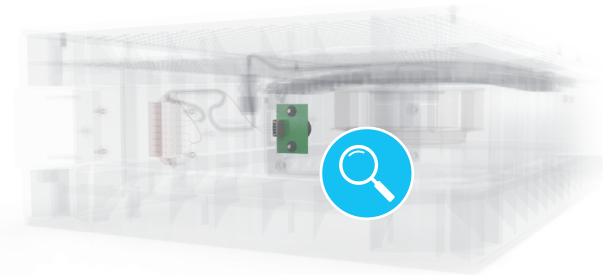

Three sensors, specifically adapted to detect of a very broad spectrum of molecules:








Formaldehyde sensor (Molecode F)
Electrochemical sensor detects formaldehyde vapours.

Acid sensor (Molecode A)
Electrochemical sensor detects inorganic acid vapours.

VOC sensor (Molecode S)
Semiconductor sensor detects volatile organic compounds.

Sensors

For optimal enclosure containment and handling safety according to your device.

 Formaldehyde (Molecode F)	 Acids (Molecode A)	 VOC (Molecode S)	Anemometer	Tachometer
Filtration performance monitoring for Formaldehyde	Filtration performance monitoring for Acids	Filtration performance monitoring for COV	Air face velocity monitoring	Monitoring of fans
Included in accordance with the eValiQuest assessment. (Automatic online analysis of chemical hazards)	Included in accordance with the eValiQuest assessment. (Automatic online analysis of chemical hazards)	Included in accordance with the eValiQuest assessment. (Automatic online analysis of chemical hazards)	 Included	 Included

FLEX® TECHNOLOGY

US patent number 7,766,732 B.

Modular filtration column

The combination of molecular and particulate filtration technologies enables one single device to be configured for the laboratory's protection needs.

This has been made possible by the design of single-dimension filter cartridges which, through vertical stacking, represent a major innovation in **Flex technology**.

The modular filtration column **adapts to the protection needs** and laboratory environment. The different Captair filtering fume hood models in the range can also be equipped with one to four filtration columns with **unrivalled** retention capacities.

This innovation from the Erlab R&D laboratory offers unprecedented **flexibility, adaptability** and **savings**. One single device can be reconfigured over time and easily reassigned to other applications.



FILTRATION CONFIGURATIONS



LIQUIDS

For dilutions, dosages, extractions, transfers, etc.



POWDERS

For sieving, grinding, weighing, formulations, compressions, etc.



LIQUIDS & POWDERS

For solutions, filtration, extractions, etc.



IN CLEAN ROOMS

For clean rooms up to ISO class 7

1C




1P



1P 1C



1C 1P

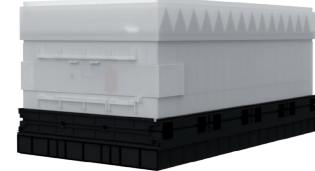

2C




2P



1P 2C



2C 1P

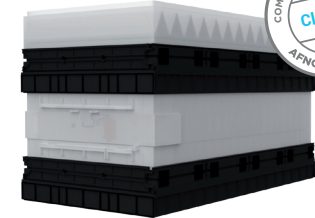

1P 1C 1P



Class 1 = Maximum safety guaranteed by the AFNOR NF X 15-211:2009 standard.
Filtration fault sensor included.



Carbon filter

Carbon Filter Module
Super-activated carbon.
Exceptional adsorption properties.



HEPA filter

HEPA Particulate Filter Module
High-efficiency filtration of particles in the air which guarantees an overall filtration efficiency of 99.995% (particles greater than 0.1 micron), according to the MPPS method in the EN 1822-1 standard.



Fan

Ventilation module
Ventilation regulation.
Low energy consumption.
Reduced noise level.

SMART TECHNOLOGY

Smart technology is a simple, innovative and intuitive method of communication for added safety.


Without the need for a screen, this communication interface allows operators, through its **simple and intuitive** operation, to focus all their attention on what is important: **handling**.

Using light and sound signals, this technology indicates the level of user protection.

Different **light and sound pulses** also inform the user of the device's operational status in real time.



THE BENEFITS OF SMART TECHNOLOGY:

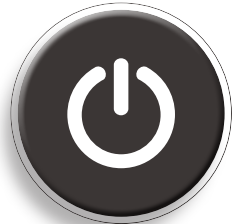

SIMPLICITY

One-touch activation for simple and intuitive operation.




SAFETY

Through the different light and sound pulses, Smart technology **keeps you informed of your protection status in real time.**



DETECTION

The exclusive detection system constantly **monitors the filtration fault status.**

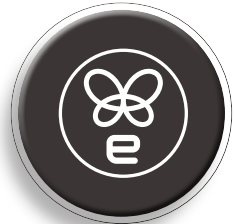

CONNECTIVITY

By connecting your device to the **eGuard** network, you can monitor the operating parameters of your Smart equipment both in real time and remotely.

1 Light pulses

Real-time communication with **flashing LED lights** provides the user with an intuitive update on the device's operational status.

2 Simplicity

One single activation key.

3 Detection system

The exclusive detection system constantly monitors filtration performance.

4 Embedded software

This service provides direct access to the following information: **the status, settings and historical data for your device.**



Smart technology alarms?

Smart technology communicates the device's functional status in real time:



Timer
(One beep sound)



Hood air face velocity
(Two beep sounds)



Prolonged opening of cabinet doors
(Two beep sounds)



Filtration fault indicator
(Three beep sounds)



Ventilation fault
(Four beep sounds)

ERLAB PRODUCTS AND THEIR APPLICATIONS



Erlab®

Flexible and connection-free filtration solutions



CHEMICAL HANDLINGS AND PROTECTION OF THE HANDLER

Captair® FILTERING CHEMICAL FUME HOODS

Captair filtering fume hoods provide unrivalled filtration quality to protect operators, and are superior to any other filtration or extraction device on the market.

By capturing toxic gases at source and filtering them through the world's most powerful filters developed by Erlab (which are 10 to 1,000 times more effective than other filters on the market), it is possible to obtain **air that is over 99.995% pure**, both inside and outside.

A system of endless benefits, providing flawless protection, not only for laboratory personnel, but also for all the company's employees, as well as everyone living nearby.

They have a powerful light communication interface, with **Smart technology** that keeps operators informed at all times about the proper functioning and level of safety of their fume hoods.

The technology offers an exceptional working comfort that removes any anxiety, allowing greater concentration on research for those who work in the fields of chemistry, pharmaceuticals, cosmetics, biochemistry, education, petrochemistry, forensics, construction, agri-food, hospitals, etc.

Flex® technology

A flexible and adaptable modular filtration column

The combination of molecular and particulate filtration technologies enables one single device to be configured for the laboratory's protection needs. This has been made possible by the design of single dimension filter cartridges which, through vertical stacking, represent a major innovation in the range. This innovation from the Erlab R&D laboratory offers unprecedented flexibility, adaptability and savings.

Sampling sound

For filter fault level detection.

(Not present on devices equipped with the Molecode S option, which automatically detects filter faults)

Electronic anemometer

This system ensures permanent control of the air face velocity, which must be between 0.4 and 0.6 m/s (in accordance with the requirements of the AFNOR NF X 15-211:2009 standard).

Easy installation

Erlab solutions are quick and easy to install. Unlike connected systems, they do not require a ventilation system linked to an air supply and extraction system to be installed. You just need one plug to get up and running, and you can get set up at any suitable moment, without complex planning.

Cable grommet


Cable grommet port.

Easy to relocate

Erlab solutions can be moved within a laboratory to suit your protection needs, without affecting the ventilation balance of the room.



CONFORMS WITH THE
AFNOR
NF X 15-211
FILTRATION STANDARD

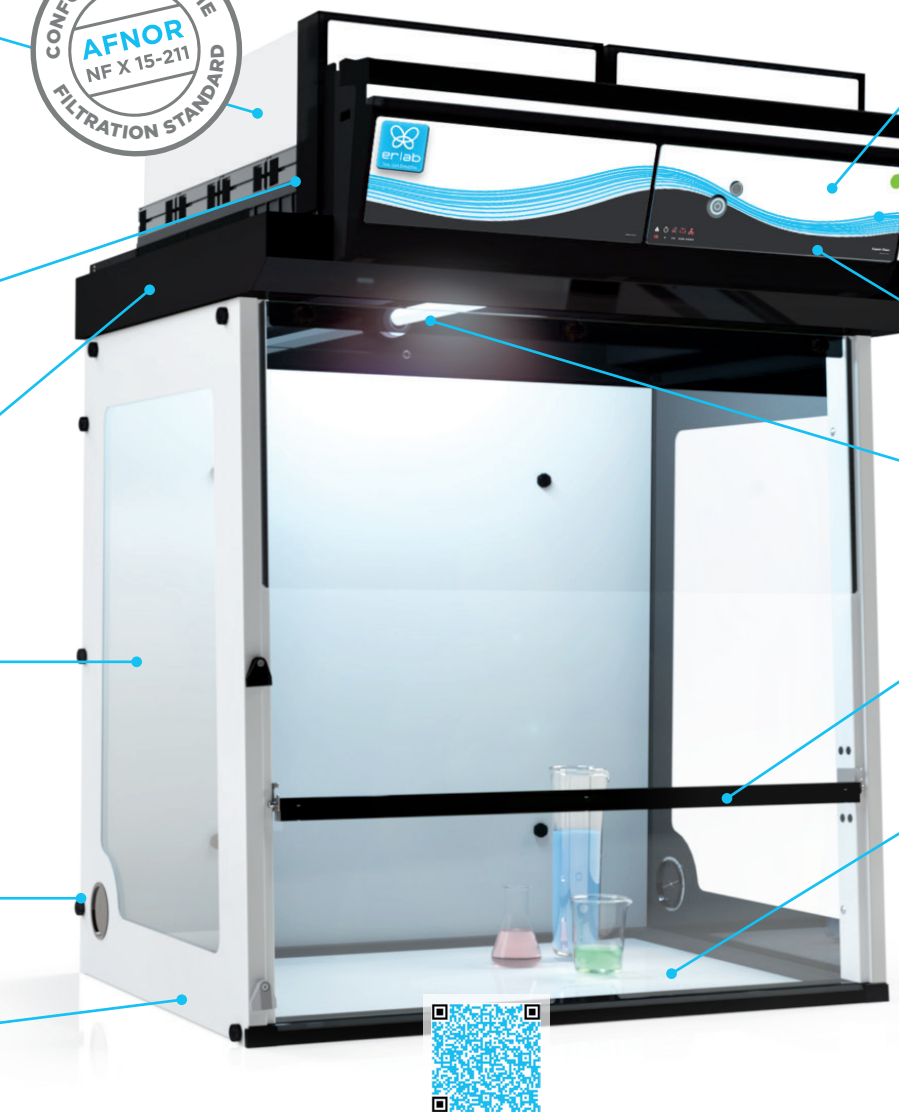

Smart technology

Light and sound pulses inform the user about their level of protection. The user will therefore be informed, in real time, of the hood's operating time, the air face velocity status, level of filtration faults and settings related to ventilation, through the various light and sound pulses.

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.

Low-energy lighting

LED lighting. One to three tubes, depending on the model. Dust and vapour proof. Bright and uniform lighting on the work surface.

New reverso panel

One panel for all working positions. With a simple movement, you can switch from a low to a medium position, in complete safety.

Interchangeable work surfaces

Glass work surface


Enamelled glass work surface with joined frame and retention tray. Ergonomic arm rest for a comfortable working position.

Phenolic resin work surface

Phenolic resin work surface with integrated retention tray and ergonomic arm rest for a comfortable working position. High chemical and mechanical resistance. Ideal for precise weighing operations.

Inox 304L work surface

Significant chemical and mechanical resistance. Rounded edges to facilitate cleaning. Integrated retention tray.

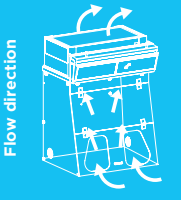

Flow direction

321

Number of filtration columns

1

Number of fans

1

Air flow rate

220 m³/h

Air face velocity

0.4 to 0.6 m/s

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

65 W

Type of opening

Oblong

Structure

Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels

Transparent and colourless PMMA with high optical purity

Filtration module

Injected polypropylene



Dimensions (mm):
Internal L 765 x D 621 x H 732*
External L 808 x D 621 x H 1110/1285**




Tested according to ASHRAE 110:1995
and conforms with BS7989
Tests and labelling C E

*Usable height of the work surface
**Max/min height of the filtration column

Possible uses



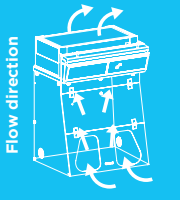

Flow direction

391

Number of filtration columns

1

Number of fans

1

Air flow rate

220 m³/h

Air face velocity

0.4 to 0.6 m/s

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

65 W

Type of opening

Oblong

Structure

Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels

Transparent and colourless PMMA with high optical purity

Filtration module

Injected polypropylene



Dimensions (mm):
Internal L 969 x D 621 x H 732*
External L 969 x D 621 x H 1110/1285**




Tested according to ASHRAE 110:1995
and conforms with BS7989
Tests and labelling C E

*Usable height of the work surface
**Max/min height of the filtration column

Possible uses



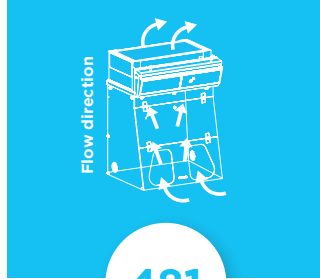

481

Number of filtration columns

1

Number of fans

1

Air flow rate
220 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
65 W

Type of opening
Oblong

Structure
Electro-galvanized and
corrosion resistant steel,
coated with thermosetting
anti-acid polymer

Side and front panels
Transparent and colourless
PMMA with high optical purity

Filtration module
Injected polypropylene




Tested according to ASHRAE 110:1995
and conforms with BS7989
Tests and labelling C E



Dimensions (mm):
Internal L 1244 x D 621 x H 732*
External L 1288 x D 621 x H 1110/1285**

*Usable height of the work surface
**Max/min height of the filtration column

Possible uses

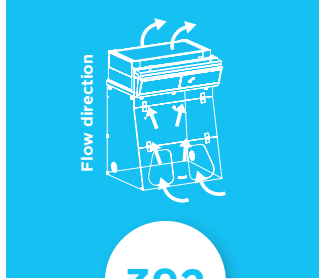

392

Number of filtration columns

2

Number of fans

2

Air flow rate
440 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
105 W

Type of opening
Reverso or sequential *see p.47*

Structure
Electro-galvanized and
corrosion resistant steel,
coated with thermosetting
anti-acid polymer

Side and front panels
Transparent and colourless
PMMA with high optical purity

Filtration module
Injected polypropylene



Tested according to ASHRAE 110:1995
and conforms with BS7989
Tests and labelling C E




Dimensions (mm):
Internal L 969 x D 749 x H 907*
External L 1005 x D 749 x H 1340/1515**

*Usable height of the work surface
**Max/min height of the filtration column

Possible uses



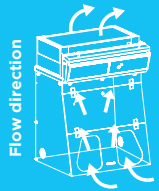

483

Number of filtration columns
3

Number of fans
3

Air flow rate
660 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
160 W

Type of opening
Reverso or sequential *see p.47*

Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Transparent and colourless PMMA with high optical purity

Filtration module
Injected polypropylene




Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 

Dimensions (mm):
Internal L 1244 x D 749 x H 907*
External L 1295 x D 960 x H 1340/1515**

*Usable height of the work surface
**Max/min height of the filtration column

Possible uses

- Liquids
- Powders
- Liquids & Powders
- Clean room

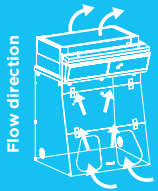

633

Number of filtration columns
3

Number of fans
3

Air flow rate
660 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
160 W

Type of opening
Reverso or sequential *see p.47*

Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Transparent and colourless PMMA with high optical purity

Filtration module
Injected polypropylene




Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 

Dimensions (mm):
Internal L 1569 x D 749 x H 907*
External L 1620 x D 960 x H 1340/1515**

*Usable height of the work surface
**Max/min height of the filtration column

Possible uses

- Liquids
- Powders
- Liquids & Powders
- Clean room

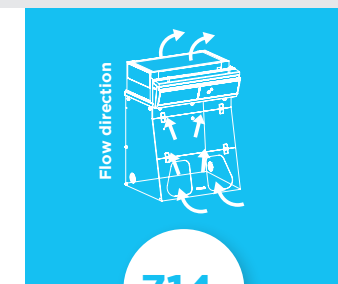

Number of filtration columns
4

Number of fans
4

Air flow rate
880 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
220 W


Type of opening
Reverso or sequential *see p.47*

Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Transparent and colourless PMMA with high optical purity

Filtration module
Injected polypropylene



Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 

Dimensions (mm):
Internal L 1769 x D 749 x H 907*
External L 1819 x D 960 x H 1340/1515**

*Usable height of the work surface
**Max/min height of the filtration column


Possible uses



EQUIPMENT

For our mobile and ductless hoods with Erlab's modular filtration columns.


Work surfaces*

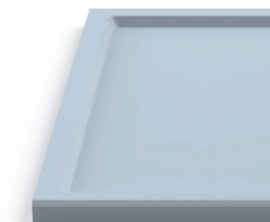

Glass work surface
Enamelled glass work surface with joined frame and retention tray.
Ergonomic arm rest for a comfortable working position.



Phenolic resin work surface
Phenolic resin work surface with integrated retention tray and an ergonomic arm rest for a comfortable working position. High chemical and mechanical resistance. Ideal for precise weighing operations.

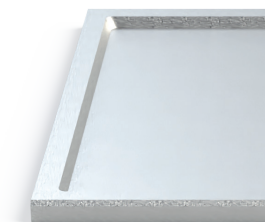

Inox 304L work surface
Significant chemical and mechanical resistance. Rounded edges to facilitate cleaning. Integrated retention tray.

*Spacing bar only for devices without a work surface

Transparent back panel



Made of synthetic glass. Offers a 360° view of handlings being carried out in the enclosure and optimises lighting.

Work bench



Mobicap
Metal work bench, equipped with four wheels, two of which are self-locking. Allows the device to be moved with complete safety.

Only available for models 321-391-392-481



Benchcap
Fixed metal work bench. Equipped with four adjustable feet allowing the unit height to be adjusted.



Benchcap optional equipment
Benchcap furniture with phenolic resin work surfaces offer multiple equipment solutions for fluids (such as a swan neck, etc.), as well as for technical gases and energy (such as nozzles on backslash, etc.) as well as an electrical socket.

Get in touch to discuss tailoring your hood to your needs.



Benchcap shelf option
Internal, semi-removeable, metal **Benchcap** shelf.

CHEMICAL HANDLINGS AND PROTECTION OF THE HANDLER

SECURE Captair® WEIGHING STATIONS

Secure Captair weighing stations guarantee users an excellent level of safety.

They offer the stability and accuracy required for applications up to 10⁻⁶g during laboratory weighing.

This is the ideal solution for high-risk operations that require a confined environment to eliminate any potential exposure of personnel to toxic products.

The Smart technology embedded in Erlab weighing stations keeps users safe: a powerful light-based communication interface keeps users informed in an intuitive way about the operating status of their weighing stations.

Flex® technology

Modular filtration technology tailored to the weighing of liquids and/or powders

The combination of molecular and particulate filtration technologies enables one single device to be configured for the laboratory's protection needs. This has been made possible by the design of single dimension filter cartridges which, through vertical stacking, represent a major innovation in the range. This innovation from the Erlab R&D laboratory offers unprecedented flexibility, adaptability and savings.

Electronic anemometer

This system ensures permanent control of the air face velocity, which must be between 0.4 and 0.6 m/s (in accordance with the requirements of the AFNOR NF X 15-211:2009 standard).

Easy to relocate

Erlab solutions can be moved easily within a laboratory to suit your protection needs, without affecting the ventilation balance of the room.

Easy installation

Erlab solutions are quick and easy to install. Unlike connected systems, they do not require a ventilation system linked to an air supply and extraction system to be installed. You just need one plug to get up and running, and you can get set up at any suitable moment, without complex planning.

Cable grommet

Cable grommet port.



Smart technology

Light and sound pulses inform the user about their level of protection. The user will therefore be informed, in real time, of the hood's operating time, the air face velocity status, level of filtration faults and settings related to ventilation, through the various light and sound pulses.

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.

Low-energy lighting

LED lighting. One to three tubes, depending on the model. Dust and vapour proof. Bright and uniform lighting on the work surface.

Waste port


It comes with a double waste bag and protective case.

Interchangeable work surfaces

Phenolic resin work surface

Phenolic resin work surface with integrated retention tray and ergonomic arm rest for a comfortable working position. High chemical and mechanical resistance. Ideal for precise weighing operations.



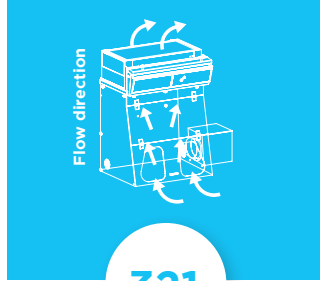

321

Number of filtration columns

1

Number of fans

1

Air flow rate

220 m³/h

Air face velocity

0.4 to 0.6 m/s

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

65 W

Type of opening

Oblong

Structure

Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer


Side and front panels

Transparent and colourless PMMA with high optical purity

Filtration module

Injected polypropylene



Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 



Dimensions (mm):

Internal L 780 x D 621 x H 732*

External L 808*** x D 621 x H 1110/1285**


*Usable height of the work surface

**Max/min height of the filtration column

*** Please note that the width does not include the waste port; see page 46

Possible uses

Liquids & Powders

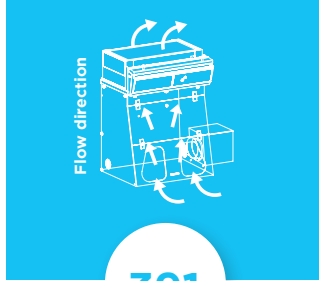

391

Number of filtration columns

1

Number of fans

1

Air flow rate

220 m³/h

Air face velocity

0.4 to 0.6 m/s

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

65 W

Type of opening

Oblong

Structure

Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels

Transparent and colourless PMMA with high optical purity

Filtration module

Injected polypropylene



Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 



Dimensions (mm):

Internal L 969 x D 621 x H 732*

External L 969*** x D 621 x H 1110/1285**

*Usable height of the work surface


**Max/min height of the filtration column

*** Please note that the width does not include the waste port; see page 46

Possible uses

Liquids & Powders



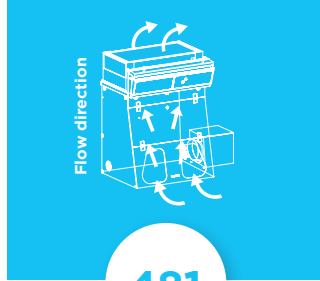

Number of filtration columns
1

Number of fans
1

Air flow rate
220 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
65 W


Type of opening
Oblong

Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Transparent and colourless PMMA with high optical purity

Filtration module
Injected polypropylene





+ Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 

Dimensions (mm):
Internal L 1244 x D 620 x H 732*
External L 1288*** x D 620 x H 1110/1285**

*Usable height of the work surface
**Max/min height of the filtration column
*** Please note that the width does not include the waste port; see page 46

Possible uses

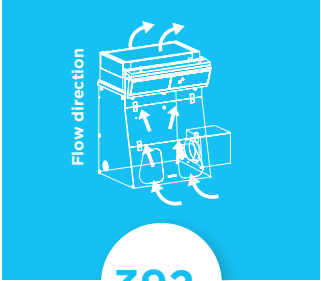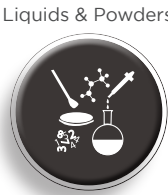

Number of filtration columns
2

Number of fans
2

Air flow rate
440 m³/h

Air face velocity
0.4 to 0.6 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
105 W


Type of opening
Reverso or sequential *see p.47*

Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Transparent and colourless PMMA with high optical purity

Filtration module
Injected polypropylene




+ Tested according to ASHRAE 110:1995 and conforms with BS7989
Tests and labelling 

Dimensions (mm):
Internal L 969 x D 749 x H 907*
External L 1005*** x D 960 x H 1340/1515**

*Usable height of the work surface
**Max/min height of the filtration column
*** Please note that the width does not include the waste port; see page 46

Possible uses

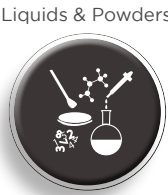

EQUIPMENT

For secure Erlab weighing stations.


Work surfaces



Phenolic resin work surface

Non-conductive material, very high mechanical and chemical resistance. Integrated retention tray. Ensures accurate weighing results and consistent operations. Prevents static charging of items introduced into the enclosure. Easy to clean.

Waste ports

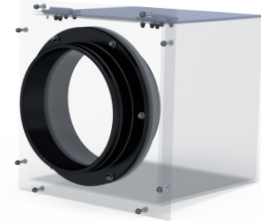

L 274 x D 278 x H 279 mm

Internal and external access secured by protective air flow. Double bag attachment system that prevents any waste from being spilled outside the enclosure. External protection against tearing or damage to the bags.

Work bench



Benchcap

Work bench that transforms the weighing unit into an independent work station. Equipped with four antivibration jacks. These allow the station height to be changed.



Benchcap optional equipment

Benchcap furniture with phenolic resin work surfaces offer multiple equipment solutions for fluids (such as a swan neck, etc.), as well as for technical gases and energy (such as nozzles on backslash, etc.) as well as an electrical socket.

Get in touch to discuss tailoring your hood to your needs.



Benchcap shelf option

Internal, semi-removeable, metal **Benchcap** shelf.

TYPE OF OPENING

For Erlab chemical handlings.


Smart 321-391-481 hoods

Oblong opening




- Front panel inclined at a 13° angle for visual comfort and a better working position.

Sequential opening

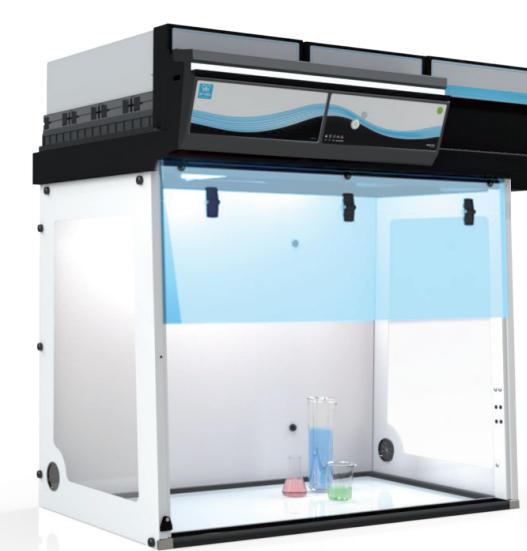

- Front panel inclined at a 4° angle for visual comfort and a better working position.

Reverso* opening

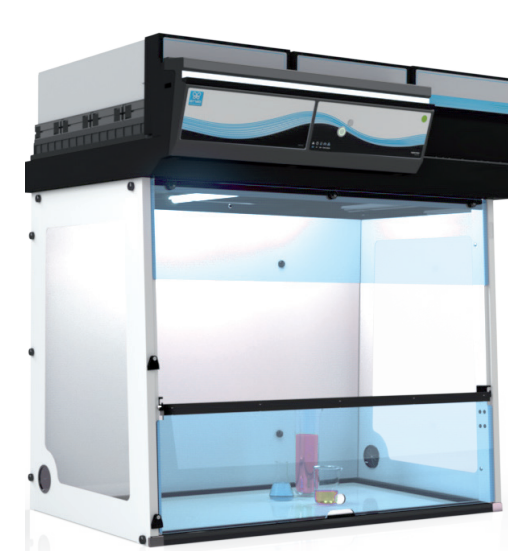

- The door can be swung open on the vertical axis for working while sitting or standing. The door can be opened completely to allow the passage of materials.
- Front panel inclined at a 4° angle for visual comfort and a better working position.

Our openings are ergonomically designed and secure to ensure that air velocity always conforms with the standard.

HANDLINGS AND SAMPLES PROTECTION

Captair® Bio PCR WORKSTATIONS

PCR workstations have a high-efficiency particulate filtration system (HEPA H14/ULPA U16) that provides a particle-free work station, eliminating all pollution around the handling. The UV lamp is also used to decontaminate the work surface and avoids biological cross-contamination between two operations.

As an option, a molecular filter can also be used to protect handlings from VOC in the laboratory air. Erlab's biological hoods are completely standalone and do not require any ventilation system connection costs nor special installation costs. With integrated **Smart technology**, the intuitive light based communication system keeps users informed about the operating status of their equipment.

Flex® technology

The combination of molecular and particulate filtration technologies enables one single device to be configured for the laboratory's protection needs. This has been made possible by the design of single dimension filter cartridges which, through vertical stacking, represent an innovation in the range. This innovation from the Erlab R&D laboratory offers unprecedented flexibility, adaptability and savings.

UV light

UV decontamination.
(Germicidal lamp - 254 nm)

Easy to relocate

Erlab solutions can be moved easily within a laboratory to suit your protection needs, without affecting the ventilation balance of the room.

Cable grommet

Cable grommet port.



Smart technology

Light and sound pulses inform the user about their level of protection. The user will therefore be informed, in real time, of the hood's operating time, the air face velocity status, level of filtration faults and settings related to ventilation, through the various light and sound pulses.

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.

Low-energy lighting

LED lighting.
Bright and uniform lighting on the work surface.

Front panels

UV-resistant front and side panels.

Work surface

Inox work surface
Significant chemical and mechanical resistance. Rounded edges to facilitate cleaning. Integrated retention tray.



320

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
25 W

Type of opening
Sequential


Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Acrylic (10 mm thick) designed to protect users from harmful UV rays and β (Beta) radiation emitted by radioactive isotopes such as: T(3H), 14C, 32P

Filtration module
Injected polypropylene




Dimensions (mm):
Internal L 765 x D 630 x H 440
External L 825 x D 630 x H 622

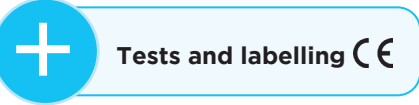

Possible uses

UV decontamination

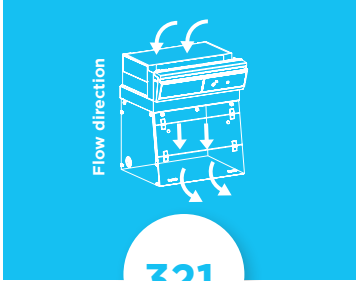

321

Number of filtration columns
1

Number of fans
1

Air flow rate*
200 / 118 m³/h of 245 / 144 m³/h

Air face velocity*
0.35 m/s

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption*
40–45 W

Type of opening
Sequential


Structure
Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels
Acrylic (10 mm thick) designed to protect users from harmful UV rays and β (Beta) radiation emitted by radioactive isotopes such as: T(3H), 14C, 32P

Filtration module
Injected polypropylene



Dimensions (mm):
Internal L 765 x D 625 x H 524
External L 810 x D 670 x H 965/1059*

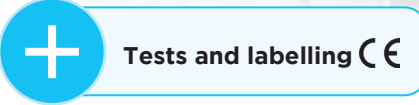

Possible uses

Ultra clean air




UV decontamination



*Min/max height depending on the filtration column

*According to the configuration of the filtration column

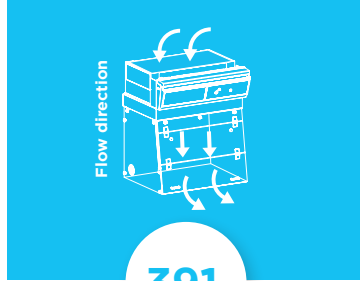

Number of filtration columns

1

Number of fans

1

Air flow rate*

200 / 118 m³/h or 245 / 144 m³/h

Air face velocity*

0.35 m/s

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption*

40–55 W

Type of opening

Sequential

Structure


Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Side and front panels

Acrylic (10 mm thick) designed to protect users from harmful UV rays and β (Beta) radiation emitted by radioactive isotopes such as: T(3H), 14C, 32P

Filtration module

Injected polypropylene

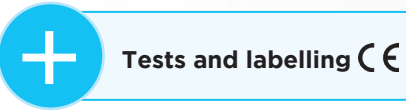

Dimensions (mm):

Internal L 969 x D 625 x H 523

External L 1013 x D 670 x H 965/1059*

*Min/max height depending on the filtration column

Possible uses

Ultra clean air




UV decontamination



EQUIPMENT

For Erlab PCR workstations.

Work surfaces

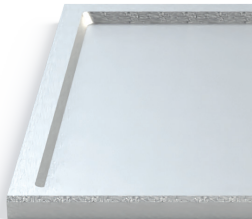

Inox 304L work surface

Significant chemical and mechanical resistance. Rounded edges to facilitate cleaning. Integrated retention tray.

Work bench



Mobicap

Metal work bench, equipped with four wheels, two of which are self-locking. Allows the device to be moved with complete safety.


Only available for models 320-321-391



Benchcap

Fixed metal work bench. Equipped with four adjustable feet allowing the unit height to be adjusted.

UV lamp

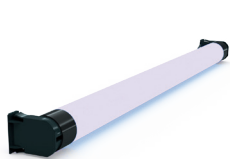

Antibacterial UV lamps:

Wave length: 254 nm

Uniform distribution of radiation on the work surface. Reflectors that eliminate shadows. Exposure adjustable from 1 to 30 minutes. Automatic UV shut-down if door is opened.



Benchcap optional equipment

Benchcap furniture with phenolic resin work surfaces offer multiple equipment solutions for fluids (such as a swan neck, etc.), as well as for technical gases and energy (such as nozzles on backsplash, etc.) as well as an electrical socket.

Get in touch to discuss tailoring your hood to your needs.



Shelf option

Internal, semi-removeable, metal **Benchcap** shelf.


Closure sensor



Sensor

The closure sensor informs the user whether the front of the device is properly closed. If the front panel is not completely closed, the UV lamp will not be triggered.

Type of opening

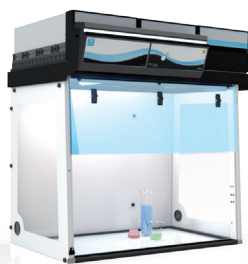

Sequential opening

The ergonomic and secure opening conforms with the standard.

HANDLINGS AND SAMPLES

PROTECTION

Captair® Flow FILTERING

CLEAN AIR WORKSTATIONS

Mobile and ductless, the filtering clean air enclosure is the ideal workstation for all handlings that require a perfectly clean environment, away from any pollution.

Erlab's high-performance filtering clean air enclosure is equipped with high-efficiency HEPA H14 (or ULPA U16) filters, which provide optimal protection against particulate contamination and also provide an ISO Class 5* work environment, according to the EN ISO 14644-1:2015 standard. Equipped with **Smart technology**, the Erlab filtering clean air enclosure allows users to focus more on their research. A powerful and intuitive light communication interface keeps them informed about their level of safety at all times. As an option, a molecular filter can also be used to protect handlings from VOC present in the laboratory atmosphere.

Flex® technology

The combination of molecular and particulate filtration technologies enables one single device to be configured for the laboratory's protection needs. This has been made possible by the design of single dimension filter cartridges which, through vertical stacking, represent a major innovation in the range. This innovation from the Erlab R&D laboratory offers unprecedented flexibility, adaptability and savings.

Low-energy lighting

LED lighting. One to three tubes, depending on the model. Bright and uniform lighting on the work surface.

Easy to relocate

Erlab solutions can be moved easily within a laboratory to suit your protection needs, without affecting the ventilation balance of the room.

Cable grommet

Cable grommet port.



Smart technology

Light and sound pulses inform the user about their level of protection. The user will therefore be informed, in real time, of the hood's operating time, the air face velocity status, level of filtration faults and settings related to ventilation, through the various light and sound pulses.

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.


Monitoring system

Continuous air flow monitoring.

Work surface

Inox work surface
Significant chemical and mechanical resistance. Rounded edges to facilitate cleaning. Integrated retention tray.



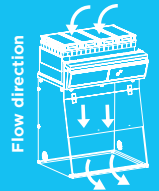

Flow direction

321

Number of filtration columns

1

Number of fans

1

Air flow rate*

320 m³/h or 150 m³/h

Air face velocity*

0.35 m/s–69 fpm

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

35–55 W

Type of opening*

Oblong or sequential

Structure

Electro-galvanised and corrosion resistant steel, coated with anti-acid polymer

Side and front panels


Chemical-resistant acrylic

Filtration module

Injected polypropylene



Dimensions (mm):
Internal L 765 x D 621 x H 732
External L 808 x D 621 x H 1110/1205**

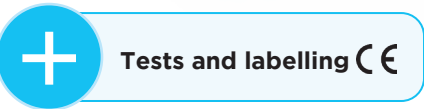


*According to the configuration of the filtration column

**Min/max height depending on the filtration column

Possible uses

Ultra clean air



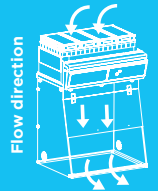

Flow direction

391

Number of filtration columns

1

Number of fans

1

Air flow rate*

345 m³/h or 150 m³/h

Air face velocity*

0.35 m/s–69 fpm

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

40–55 W

Type of opening*

Oblong or sequential

Structure

Electro-galvanised and corrosion resistant steel, coated with anti-acid polymer

Side and front panels


Chemical-resistant acrylic

Filtration module

Injected polypropylene



Dimensions (mm):
Internal L 969 x D 621 x H 732
External L 1013 x D 621 x H 1110/1205**

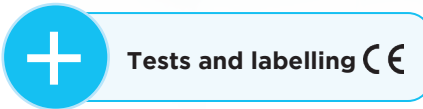


*According to the configuration of the filtration column

**Min/max height depending on the filtration column

Possible uses

Ultra clean air



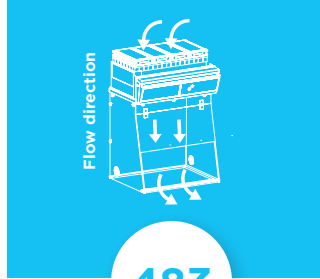

Number of filtration columns
3

Number of fans
3

Air flow rate*
770 m³/h or 530 m³/h

Air face velocity*
0.35 m/s–69 fpm

Voltage/frequency
110–230 V / 50–60 Hz


Power consumption
90–95 W

Type of opening
Sequential

Structure
Electro-galvanised and corrosion resistant steel, coated with anti-acid polymer

Side and front panels
Chemical-resistant acrylic


Filtration module
Injected polypropylene

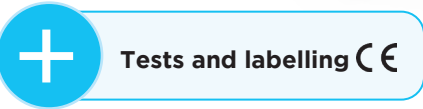

Dimensions (mm):
Internal L 1244 x D 749 x H 907
External L 1295 x D 960 x H 1340/1435**

Possible uses

Ultra clean air

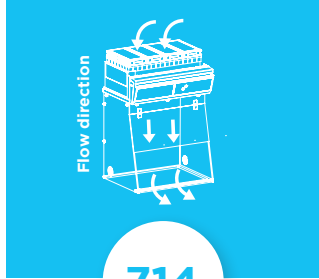

Number of filtration columns
4

Number of fans
4

Air flow rate*
1040 m³/h or 690 m³/h

Air face velocity*
0.35 m/s–69 fpm

Voltage/frequency
110–230 V / 50–60 Hz




Power consumption
110 W

Type of opening
Sequential

Structure
Electro-galvanised and corrosion resistant steel, coated with anti-acid polymer

Side and front panels
Chemical-resistant acrylic

Filtration module
Injected polypropylene

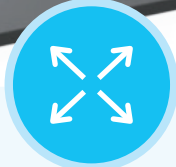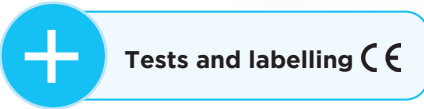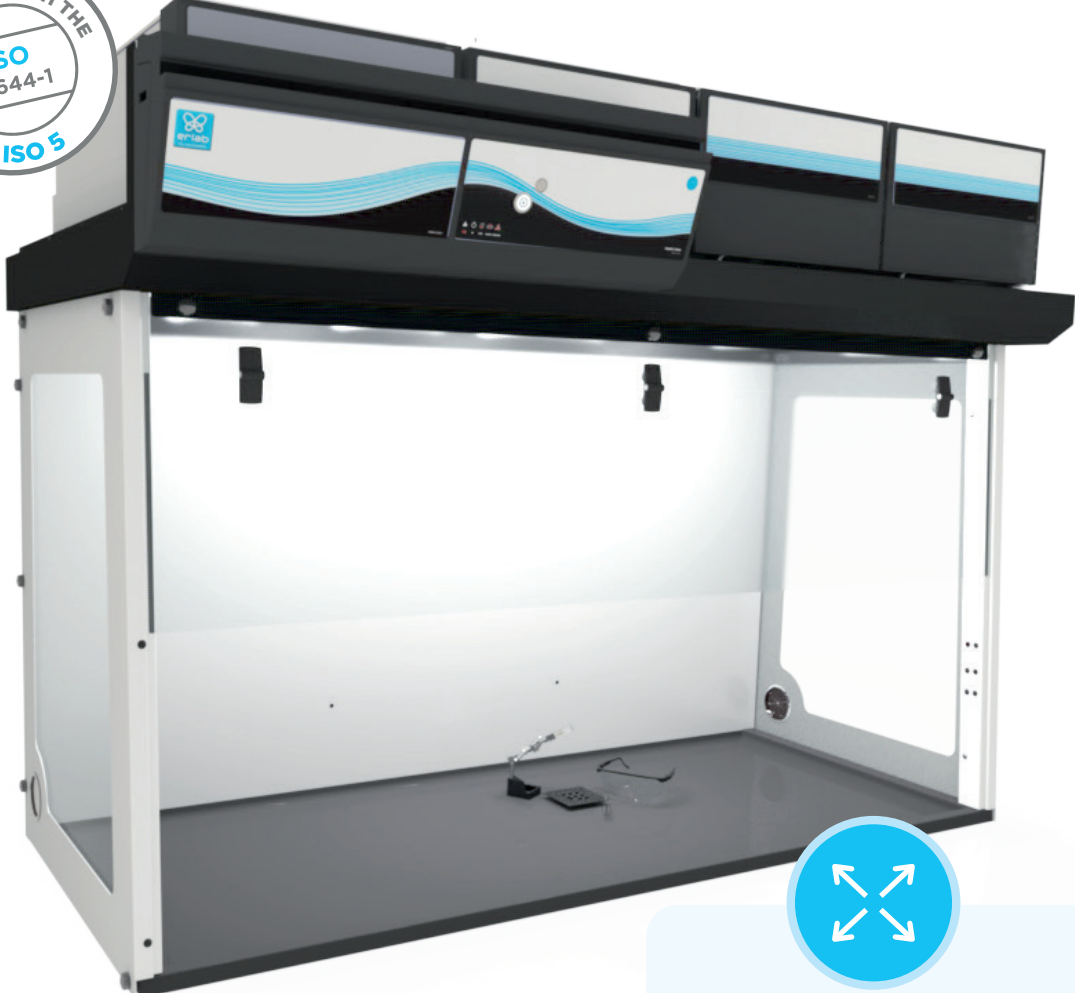

Dimensions (mm):
Internal L 1769 x D 749 x H 902
External L 1819 x D 960 x H 1340/1435**

Possible uses

Ultra clean air



*According to the configuration of the filtration column

**Min/max height depending on the filtration column

*According to the configuration of the filtration column

**Min/max height depending on the filtration column

EQUIPMENT


For Erlabfiltering clean air workstations.

Work surfaces



Phenolic resin work surface

Phenolic resin work surface with integrated retention tray and an ergonomic arm rest for a comfortable working position. High chemical and mechanical resistance. Ideal for precise weighing operations.

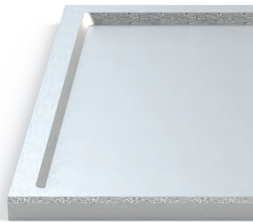

Inox 304L work surface

Significant chemical and mechanical resistance. Rounded edges to facilitate cleaning. Integrated retention tray.

Work bench



Mobicap

Metal work bench, equipped with four wheels, two of which are self-locking. Allows the device to be moved with complete safety.

Only available for models 321-391




Benchcap optional equipment

Benchcap furniture with phenolic resin work surfaces offer multiple equipment solutions for fluids (such as a swan neck, etc.), as well as for technical gases and energy (such as nozzles on backsplash, etc.) as well as an electrical socket.

Get in touch to discuss tailoring your hood to your needs.


Type of opening

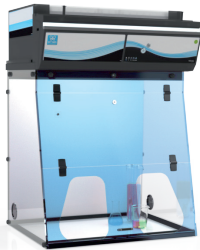

Smart Flow 321-391 fume hoods

Oblong opening

The ergonomic and secure opening conforms with the standard.

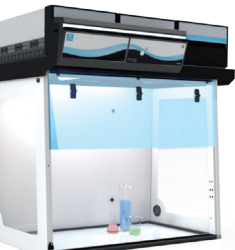



Smart Flow 483-714 fume hoods

Sequential opening

The ergonomic and secure opening conforms with the standard.

WORK SURFACES

For Erlab fume hoods.

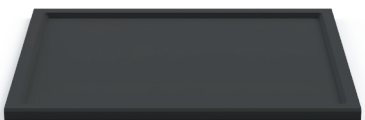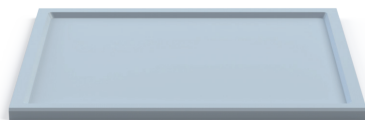

Work surfaces

	Glass work surface	Inox 304L work surface	Phenolic resin work surface
Solvents	+++	+++	+++
Acids	+++	-	++
Excluding hydrofluoric acid	×	×	-
Sodium hypochlorite (bleach)	+++	×	+
Impact resistant	+	++	+++
Scratch resistant	+++	-	-
Load capacity	50 kg/m²	110 kg/m²	110 kg/m²
Load resistant	-	+	++
Ease of cleaning	++	+++	+++
Stability	+	-	++
Equipment available	N/A	N/A	Sinks, taps (contact us for more information)

CHEMICAL STORAGE

Captair® FILTERING STORAGE CABINETS

Filtering storage cabinets provide excellent filtration quality to protect users from chemical hazards. Their technology is suitable for all types of chemicals.

Equipped with transparent doors, they offer users an optimal view of the products stored inside, and warn them in case of prolonged opening to limit inhalation risks. Erlab's new range of connected cabinets employ **Smart technology** as a simple and innovative means of communication. Using light pulses and being simple and intuitive to use, this communication interface allows operators to focus all their attention on what is important: their work.

Flex® technology

The combination of molecular and particulate filtration technologies enables one single device to be configured for the laboratory's protection needs. This has been made possible by the design of single dimension filter cartridges which, through vertical stacking, represent a major innovation in the range. This innovation from the Erlab R&D laboratory offers unprecedented flexibility, adaptability and savings.

Easy to relocate

Erlab solutions can be moved easily within a laboratory to suit your protection needs, without affecting the ventilation balance of the room.

Safe storage

Store your chemicals securely with a lock and key system.

Low consumption

Very low energy consumption limits the operational cost.

Storage

Transparent doors offer an optimal view of the storage space. Multiple ergonomic storage solutions, including removable shelves, pull-out doors and drawers.



Smart technology

Light and sound pulses inform the user about their level of protection. Different light and sound pulses also inform them of the device's operational status in real time.

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard


Connect to the eGuard interface to control and monitor your protection settings remotely.

Detection

The exclusive detection system constantly monitors the filtration fault status.

Door sensor

For optimal safety, opening the door triggers an increase in air velocity. In cases of prolonged opening, the alarm is triggered to warn the user and limit their risk of inhalation.

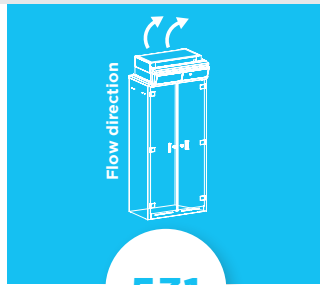

531

Number of fans
1

Air flow rate
10 m³/h

Voltage/frequency
110-230 V / 50-60 Hz

Power consumption
12 W

Doors
Transparent and colourless PMMA with high optical purity. Strong resistance to numerous aggressive chemical agents


Storage capacity
Sliding doors:
About 72 x 20 ml bottles or
About 10 x 500 ml bottles

Structure
Corrosion resistant and electro-galvanized steel, coated with thermosetting anti-acid polymer


Filtration module
Injected polypropylene



Dimensions (mm):
External L 506 x P 242 x H 406

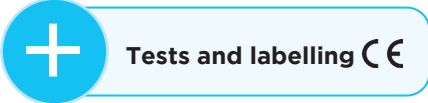

Possible uses

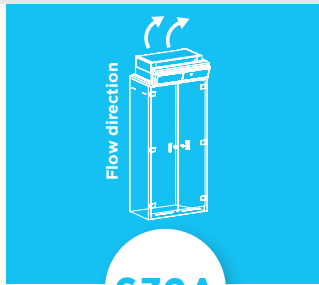

632A

Number of fans
1

Air flow rate
10 m³/h

Voltage/frequency
110-230 V / 50-60 Hz

Power consumption
16 W

Doors
Transparent and colourless PMMA with high optical purity. Strong resistance to numerous aggressive chemical agents


Storage capacity
Sliding doors:
About 24 x 500 ml bottles (high shelves) / About 20 x 1L bottles (low shelves)

Structure
Corrosion resistant and electro-galvanized steel, coated with thermosetting anti-acid polymer


Filtration module
Injected polypropylene



Dimensions (mm):
External L 603 x P 471,2 x H 685

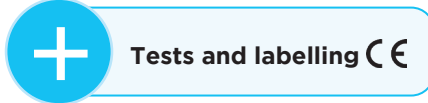

Possible uses

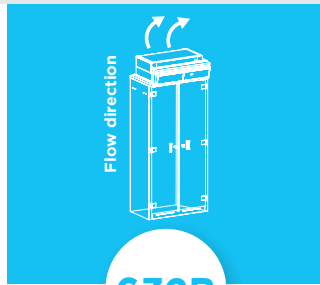

632B

Number of fans
1

Air flow rate
11 m³/h

Voltage/frequency
110-230 V / 50-60 Hz

Power consumption
20 W

Doors
Transparent and colourless
PMMA with high optical purity.
Strong resistance to numerous
aggressive chemical agents


Storage capacity
Sliding doors:
About 44 x 1L bottles
Sliding shelf :
About 25 x 1L bottles de 1L (tray)
About 35 x 100 ml bottles (shelf)
About 70 x 50 ml bottles (shelf)

Structure
Corrosion resistant
and electro-galvanized steel,
coated with thermosetting
anti-acid polymer


Filtration module
Injected polypropylene



Dimensions (mm):
External L 654 x P 583,6 x H 715

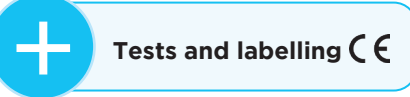

Possible uses
Liquids & Powders

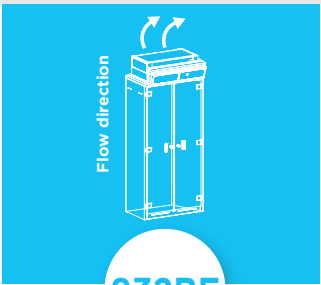

832BF

Number of fans
1

Air flow rate
70 m³/h

Voltage/frequency
110-230 V / 50-60 Hz

Power consumption
20 W

Doors
Transparent and colourless
PMMA with high optical purity.
Strong resistance to numerous
aggressive chemical agents.

Storage capacity
Double doors/shelves:
120 glass bottles of 1L

Structure
Corrosion resistant
and electro-galvanized steel,
coated with thermosetting
anti-acid polymer

Filtration module
Injected polypropylene



Dimensions (mm):
External L 900 x D 500 x H 2060




Possible uses
Liquids



Powders



Possible uses

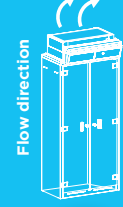

Flow direction

834

Number of fans

1

Air flow rate

220 m³/h

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

45 W

Doors

Transparent and colourless PMMA with high optical purity. Strong resistance to numerous aggressive chemical agents.

Storage capacity

Pull-out doors:
100 bottles of 1L

Double doors/shelves:
120 bottles of 1L

Structure

Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Filtration module

Injected polypropylene



FILTER PERFORMANCE
AFNOR
NF X 15-211
TESTED TO COMPLY WITH






Dimensions (mm):
External L 900 x D 660 x H 2140




Tests and labelling €€

Possible uses


Liquids




Powders



Liquids & Powders



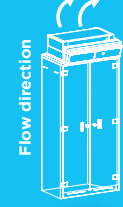

Flow direction

1634

Number of fans

1

Air flow rate

220 m³/h

Voltage/frequency

110–230 V / 50–60 Hz

Power consumption

45 W

Doors

Transparent and colourless PMMA with high optical purity. Strong resistance to numerous aggressive chemical agents.

Storage capacity

Pull-out doors:
200 glass bottles of 1L

Double doors/shelves:
240 glass bottles of 1L

Structure

Electro-galvanized and corrosion resistant steel, coated with thermosetting anti-acid polymer

Filtration module

Injected polypropylene



FILTER PERFORMANCE
AFNOR
NF X 15-211
TESTED TO COMPLY WITH






Dimensions (mm):
External L 1805 x D 660 x H 2175/2270*




Tests and labelling €€

Possible uses


Liquids



Powders



Liquids & Powders



*Min/max height depending on the filtration column

EQUIPMENT

For Erlab filtering storage cabinets.

Shelves




Shelves

Fixed and height-adjustable shelves (depending on the height of the bottles)

Integrated retention tray. Four-litre capacity.

Dimensions: 360 x 380 mm


HCL sealed box

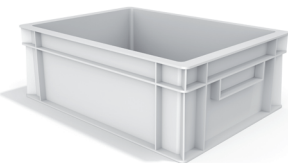

HCL sealed box

Sealed box for storing strong acids.
(store box at the bottom of the cupboard)

Keys

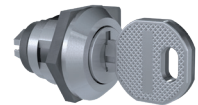

Secure storage with a lock and key system

Store your chemicals securely
with a lock and key system.

Absorbent mat



Absorbent mat

In addition to the retention tray on the
shelf, the mat holds various chemicals to
prevent them from spreading.

TYPE OF OPENING

For Erlab filtering storage cabinets.

Smart cabinets

Double doors



Pull-out doors



CHEMICAL STORAGE

Chemtrap® FILTRATION BOXES FOR FIRE-RESISTANT CABINETS

Filtration boxes allow laboratories with a safety cabinet to combine fire protection with inhalation protection. Now compatible with many safety cabinets (Asecos, Duperthal, etc.) that have an air extraction diameter of 75, 80 or 100 mm. They eliminate the harmful fumes emitted from the products stored inside the cabinet*. Thanks to their filtered air recirculation principle, they contribute to the purification of the ambient air. They require no ductwork and can be connected quickly and easily by connecting to your safety cabinet.

There are two filtration boxes available:

H402: for fire-resistant safety cabinets and high or low storage cabinets.

V201: for fire-resistant safety cabinets and storage under a bench.

Filtration

The filtration box contributes to the purity of the air in the laboratory. Filtering the air inside the storage cabinet, it captures toxic vapours at source and circulates clean air through the laboratory, contributing to a better and safer working environment.

Safety cabinet

The system transforms any fire-resistant cabinet or other chemical storage cabinet into a safety cabinet with standalone filtration. No ductwork is required.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.

Easy to relocate

Quick and easy installation.
Small space requirement. **(V201)**

Quick and easy installation. Handle for easier transport and installation. **(H402)**



H402



V201

Smart technology

Light and sound pulses inform the user about their level of protection. Different light and sound pulses also inform them of the device's operational status in real time.

Detection

The exclusive detection system constantly monitors the filtration fault status.

Connection

A flexible air duct connects the filtration box and the fire-resistant cabinet. **(V201)**


Four connection options: left, right, at the back or underneath. **(H402)**

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

Low consumption

Very low energy consumption.
Reduces running costs.

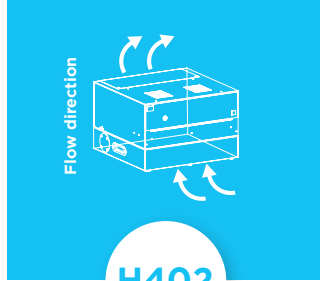

H402

Number of fans
1

Air flow rate
60 m³/h

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
20 W

Flexible duct connection
1 metre

Structure
Coated and corrosion resistant
steel (100% polyester)

Filtration module
Injected polypropylene




Dimensions (mm):
Open L 490/542* x D 520 x H 464
Closed L 490/542* x D 520 x H 326

*Min/max width depending on the type of connection



Tests and labelling C €

Possible uses

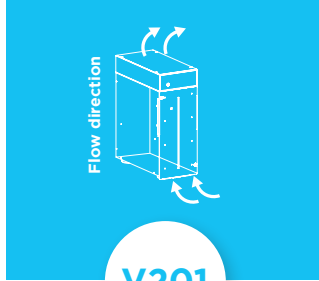

V201

Number of fans
1

Air flow rate
35 m³/h

Voltage/frequency
110–230 V / 50–60 Hz

Power consumption
15 W

Flexible duct connection
1 metre

Structure
Coated and corrosion resistant
steel (100% polyester)

Filtration module
Injected polypropylene



Dimensions (mm):
Open L 245 x D 1004 x H 632
Closed L 245 x D 554 x H 632



Tests and labelling C €

Possible uses



FILTRATION BOXES

There are two versions of Erlab filtration boxes available.

H402

For high or low cabinets




Connection: at the back, right or left.

EQUIPMENT

For Erlab filtration boxes

Flexible duct connection

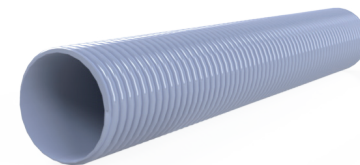

Filtration box duct connection

Adjustable and flexible air duct that connects easily to your storage cabinet's air outlet.

FILTRATION BOXES

There are two versions of Erlab filtration boxes available.

V201

For cabinets under a bench




Connection: at the back.

EQUIPMENT

For Erlab filtration boxes

Flexible duct connection

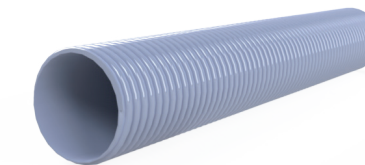

Filtration box duct connection

Adjustable and flexible air duct that connects easily to your storage cabinet's air outlet.

HALO AIR PURIFIERS

(Ceiling installation)

Halo AIR PURIFIERS

Halo air purifiers, installed on the ceiling, capture pollutants in the laboratory air or other rooms at source and trap them permanently in molecular filters or high efficiency particulate filters. An ideal solution for reducing air renewal costs, Erlab's air purifiers ensure excellent air quality in laboratories or other rooms of your building, without the need for unwarranted air renewal and the ensuing high levels of energy consumption. Available in four versions so as to filter many pollutants, Erlab's **Halo** air purifiers monitor the air quality of laboratories or other rooms where it is installed, in real time and filter the ambient air, in order to maintain an optimal level of hygiene and safety.

As a standalone system, completely independent of a building's ventilation system, Erlab's air purifiers can be easily installed, in the ceiling, in any buildings. Equipped with **Smart technology**, they allow you to monitor your room's air quality remotely, through your mobile phone. This connectivity also means that you benefit from exclusive guarantees and services.

Detection

A dedicated air quality sensor (HALO 35 only).

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.

Clean air diffusers

Four clean air diffusers provide uniform circulation of air throughout the room.

Air quality

Air purifiers provide a high level of air quality without having to connect to a central HVAC system. They also generate substantial energy savings.

Low consumption

Very low energy consumption.
Reduces running costs.



HALO 35



Smart technology

Light and sound pulses inform the user about their level of protection. Different light and sound pulses also inform the user about the ambient air quality in real time.

Ceiling installation


Safer and more effective. Both good reasons for placing the air purification station in the center of your ceiling. In this location it will also be protected from the movement of people in the room and situated among the airborne droplets that mainly circulate in the upper part of the room.

Easy access

Easy access to the chamber to replace the filters.

Filtration

Air purifiers guarantee the decontamination of the air in any laboratory, office or classroom against all chemical, viral and bacterial pollutants.

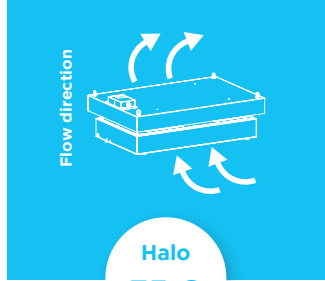

- Number of fans**
1
- Air flow rate**
220 m³/h
- Voltage/frequency**
110–230 V / 50–60 Hz
- Power consumption**
50 W
- Mode of operation**
24/7, Day/Night
Min/max detection,
detection value only
- Weight (kg)**
31 kg (including filter)
- Volume protected**
55 m³ or a surface area
of approximately 25 m² with
a ceiling height of 2.20m
- Structure**
Electro-galvanized and
corrosion resistant steel, coated
with thermosetting
anti-acid polymer
- Filtration module**
Injected polypropylene





+ Tests and labelling C €



Dimensions (mm):
External L 592 x D 892 x H 303

Possible uses

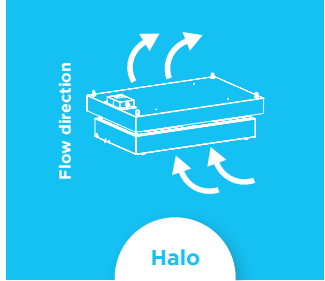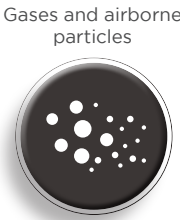

- Number of fans**
1
- Air flow rate**
300 m³/h
- Voltage/frequency**
110–230 V / 50–60 Hz
- Power consumption**
50 W
- Mode of operation**
24/7, Day/Night
- Weight (kg)**
31 kg (including filter)
- Volume protected**
70 m³ or a surface area
of approximately 30 m² with
a ceiling height of 2.50m
- Structure**
Electro-galvanized and
corrosion resistant steel,
coated with thermosetting
anti-acid polymer
- Filtration module**
Injected polypropylene




+ Tests and labelling C €



Dimensions (mm):
External L 592 x D 892 x H 303

Possible uses

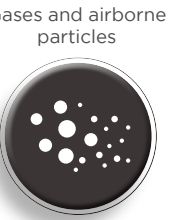

SMART TECHNOLOGY

Simple intuitive communication by light ring pulsations shows the status of the HALO unit


Alarm contact for BMS & BAS systems

3 molecular filters possible

Dedicated air quality sensor

Easy filter & fan replacement access chamber

Prefilter

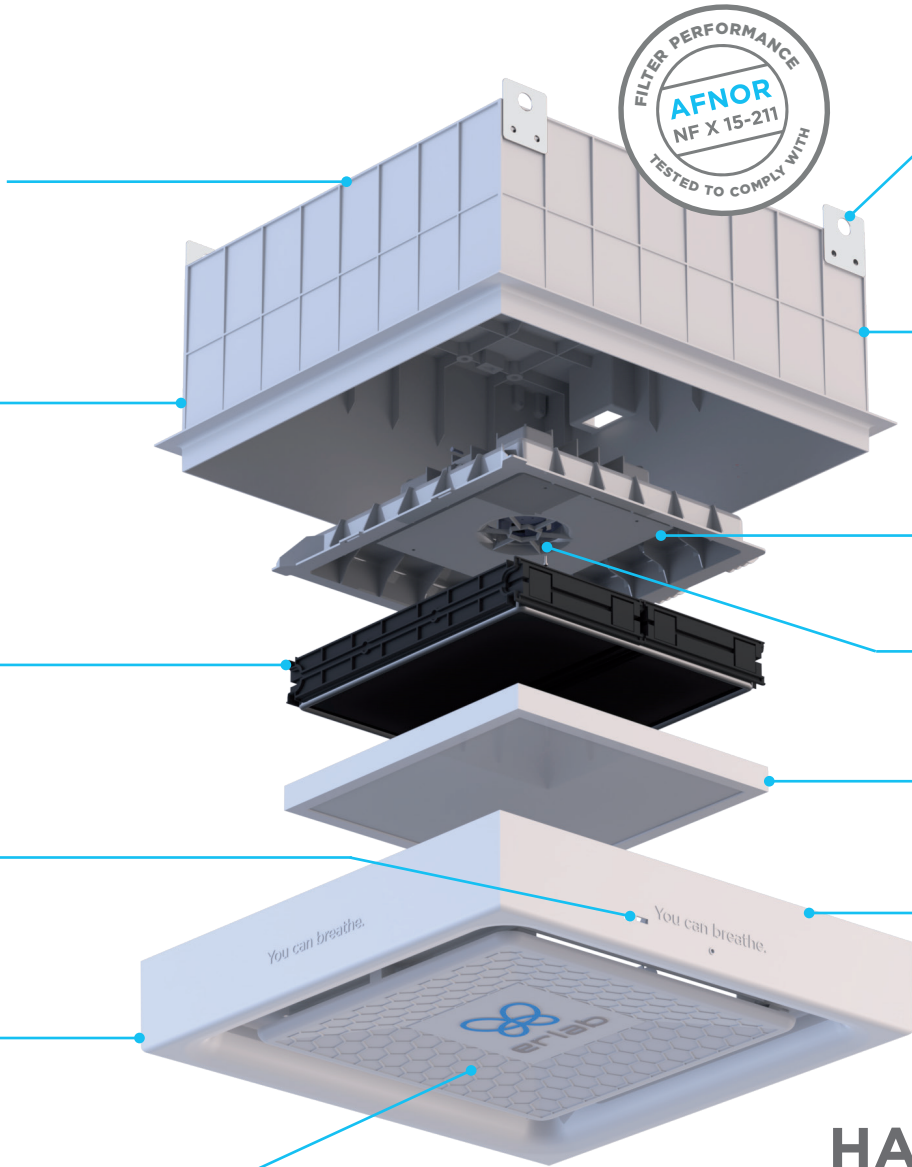

2 types of fixing :
4 plates or 4 brackets

Ethernet port for remote safety monitoring

Particulate post-filter (Molecular or Double-filtration)

Pulse-Width Modulation controlled fan

Particulate filter H14

Clean air injectors diffuse filtered air evenly back into the room



HALO 25



Halo
25 C

Number of fans
1

Air flow rate
110 m³/h

Voltage/frequency
110-230 V 50/60 Hz

Power consumption
20 W

Mode of operation
24/24h - 7/7, Night/Day, Min Max detection, Detection value only

Weight (kg)
17,5 kg (filter included)

Volume protected
22,5m³ or a surface area of 9m² with a ceiling height of 2m50

Structure
ABS (Acrylonitrile Butadiene Styrene) / Injected polypropylene

Filtration module
Injected polypropylene



Dimensions (mm):
External L 615 x P 615 x H 350



Tests and labelling C €

Gases and airborne particles



Possible uses



Halo
25 P

Number of fans
1

Air flow rate
150 m³/h

Voltage/frequency
110-230 V 50/60 Hz

Power consumption
20 W

Mode of operation
24/24h - 7/7, Night/Day, Min Max
detection, Detection value only

Weight (kg)
14,5 kg (filter included)

Volume protected
30m³ or a surface area of 12m²
with a ceiling height of 2m50

Structure
ABS (Acrylonitrile Butadiene
Styrene) / Injected polypropylene

Filtration module
Aluminum



Dimensions (mm):
External L 615 x P 615 x H 350

+ Tests and labelling C €

Possible uses

Gases and airborne
particles



Halo 25
Bifiltration

Number of fans
1

Air flow rate
110 m³/h

Voltage/frequency
110-230 V 50/60 Hz

Power consumption
35 W

Mode of operation
24/24h - 7/7, Night/Day, Min Max
detection, Detection value only

Weight (kg)
19 kg (filter included)

Volume protected
22,5m³ or a surface area of 9m²
with a ceiling height of 2m50

Structure
ABS (Acrylonitrile Butadiene
Styrene) / Injected polypropylene

Filtration module
Injected polypropylene /
Aluminum



Dimensions (mm):
External L 615 x P 615 x H 350

+ Tests and labelling C €

Possible uses

Gases and airborne
particles



REMOVAL OF AIR POLLUTION

Halo Sense AIR QUALITY SENSOR

Halo Sense air quality sensors monitor the air quality in your laboratory in real time. Available in **three different versions**, they are able to detect a wide spectrum of pollutants and vapours, including Volatile Organic Compounds (**VOC**), **acid** vapours and **Formaldehyde** vapours. If pollution is detected, the sensors alert the user immediately via light and sound pulses.

Halo Sense air quality sensors are very easy to install in all laboratories. Equipped with **Smart technology**, they allow you to monitor your laboratory's air quality remotely, through your mobile phone. This connectivity also means that you benefit from exclusive guarantees and services.

Detection

There are two Halo Smart Sense versions:

- VOC (Volatile Organic Compounds)
- Acids
- Formaldehyde

Air quality

The Halo Sense Smart continuously monitors the air quality in your laboratory and detects ambient air pollution for a wide spectrum of pollutants and vapours.

Easy to relocate

Quick and easy installation.
Installation via removable turntable: wall, bench, ceiling

Energy savings

The Halo Sense Smart can communicate with building management systems. If there is no pollution, it is possible to reduce overconsumption of clean air and the resulting energy consumption.



Smart technology


Light and sound pulses inform the user about their level of protection. These light and sound pulses also inform the user about the air quality in the laboratory in real time.

Embedded software

Connect your device directly to your computer using an RJ45 cable to access information on your device.

eGuard

Connect to the eGuard interface to control and monitor your protection settings remotely.

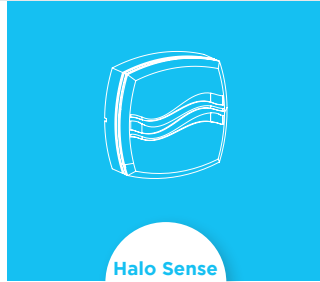

Halo Sense
VOC

Case
ABS

Power consumption
10 W

Voltage/frequency
110–230 V / 50–60 Hz

Detection
Semiconductor sensor

Weight (kg)
1 kg


Sensitivity
Five settings

Coverage
between 14 and 40 m²


Structure
Acrylonitrile Butadiene
Styrene (ABS)



Dimensions (mm):
External L 181 x W 181 x H 67

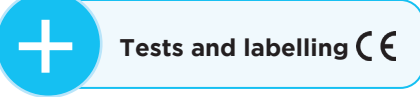

Possible detections

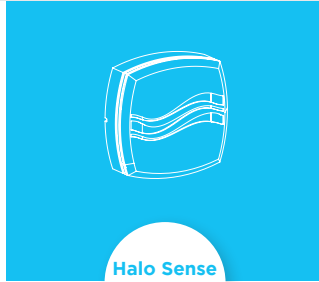

Halo Sense
Acids

Case
ABS

Power consumption
10 W

Voltage/frequency
110–230 V / 50–60 Hz

Detection
Electrochemical sensor

Weight (kg)
1 kg


Sensitivity
Three settings

Coverage
between 14 and 40 m²

Structure
Acrylonitrile Butadiene
Styrene (ABS)




Dimensions (mm):
External L 181 x W 181 x H 67

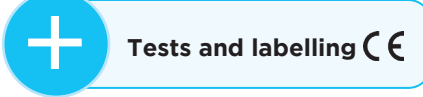

Possible detections



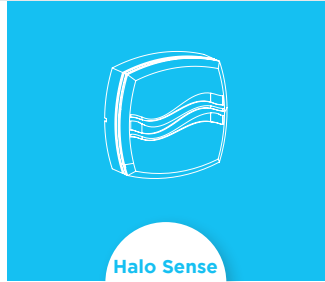

Halo Sense
Formaldehyde

Case
ABS

Power consumption
10 W

Voltage/frequency
110–230 V / 50–60 Hz

Detection
Electro-chemical sensor

Weight (kg)
1 kg

Sensitivity
Five settings

Coverage
between 14 and 40 m²

Structure
Acrylonitrile Butadiene
Styrene (ABS)



Dimensions (mm):
External L 181 x W 181 x H 67



Tests and labelling C E

Possible detections

Formaldehyde
For the detection of
formaldehyde vapours



Halo Sense monitors the air quality
in your laboratory in real time.



VOC

For the detection of
Volatile Organic Compounds



Formaldehyde

For the detection of
formaldehyde vapours



Acids

For the detection
of acid vapours

PROTECTION OF INVESTIGATIONS

TRANSPORTABLE AND DISPOSABLE Pyramid® GLOVE ENCLOSURE

Made of high-quality transparent PVC, the **Pyramid isolation enclosure** is a multi-functional, transportable and disposable "glove bag" for personal protection, as well as for the investigation and transportation of biological elements. Easily transportable, lightweight and assembled in a few seconds, it allows users to intervene anywhere, while being fully protected.

The inclined structure provides an ergonomic working position. Fully waterproof, it is ideal for any kind of field investigation, as well as for highly secure use in the laboratory. Each enclosure is highly resilient and is tested before being put on the market. Compressed air inflation is performed at a pressure of 2.5 mm Hg. A compliance certificate is provided. Each enclosure can be identified by its serial number.

Safety

The Erlab glove bag is made of high-quality, resistant and transparent PVC and is designed for the protection of people, as well as the safe investigation and transportation of biological elements. To ensure maximum safety, each enclosure is tested before being put on the market. Compressed air inflation is performed at a pressure of 2.5 mm Hg. A compliance certificate is provided with each enclosure, which is identified by a serial number.

Performance

Its quality construction makes it a vital tool for handling potentially contaminated objects, taking samples to different sites, working in an inert atmosphere, revealing fingerprints and, of course, protecting people from splashes during biopsies. Fully watertight and resistant, it keeps handlings free of dust and humidity.

Gloves

PVC medical gloves, assembled by high-frequency welding

Easy transportation

High-resistance grommet for transporting and suspending the enclosure.

Optimal impermeability

Zip opening and two-groove closure for an optimal seal.


Filling valve

Valve for potential filling of the enclosure with an inert gas (Nitrogen). Inflation of the enclosure for normal use up to 2.5 mm Hg.

Easy to install and use

Assembled in seconds, can be used anywhere. The inclined structure provides an ergonomic working position. Transportable and light, this safety enclosure is suitable for many laboratory and field requirements. It is equipped with a flexible PVC packaging for easy transportation.



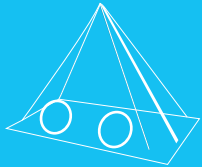

Pyramid

Enclosure and base

Flexible PVC, assembled by high-frequency welding

Closure

Double-sealed groove
(600 mm)

Medical gloves

In butyl and PVC sleeves

Usage temperature

-25°C to +45°C

Weight (kg)

1.55 kg



Dimensions (mm):

External L 860 x D 560 x H 725

Possible applications:

- Handling potentially contaminated objects
- Taking samples to different sites
- Working in an inert atmosphere
- Revealing fingerprints
- Protecting people from splashes during biopsies
- Keeping handlings free of dust and humidity

USAGE EXAMPLES

For the Erlab glove bag.



Handling on
a car



Handling on
a wall



Handling on
a work surface



Storage of several enclosures
on a crossbar before analysing
samples in the laboratory

SERVICES



eGuard allows you to monitor all of your devices remotely, and increases the safety of all users.

You can use the eGuard interface to connect your devices, including fume hoods, GreenFumeHood 3 fume cupboards, filtering storage cabinets and any other Erlab-connectable device. In just one glance, you can control the proper functioning of your equipment and the air quality in your laboratory remotely. That's our guarantee of enhanced safety.

Discover the following functions:

- Consulting usage data
- Historical data access
- Alerts
- Control over all your devices

Statistics you can access:

- **Filtration:** displays the number of events related to filtration quality.
- **Door/front opening:** displays the door or front opening percentage and the number of times the device is in containment alarm mode.
- **Temperature:** trend graph showing the temperature in the enclosure.
- **Humidity:** trend graph showing the humidity in the enclosure.



The selection of protective equipment depends on the risk to which the operator is exposed. This selection cannot be made without an in-depth risk assessment.

The **eValiQuest** service supports you throughout this process.

An expert chemist will help you to fill out the initial questionnaire to specify the handlings that you plan to carry out. Within 48 hours, our laboratory specialists will inform you of the type of device and filtration technology suitable for your use. **We are committed to ensuring your protection by certifying the feasibility of your handlings.**

A **ValiPass** usage certificate will indicate the specific chemicals used, the type of filter and the estimated service life for which your device has been validated.

What is the eValiQuest service?

This digital tool for assessing the chemical hazards associated with your handlings enables you to configure the most suitable filtering fume hood for your activity. Our chemical experts will recommend the best solution which combines optimal equipment service life and your protection.



NEW

Discover our new eValiQuest service on

www.evaliquest.erlab.com



ERLAB MAINTENANCE

Installation, Training and Maintenance Services.

Simplified maintenance management

Assisted by our network of maintenance partners, Erlab Maintenance guarantees you optimal monitoring of all your devices.

This is so you can focus on your activities in the best possible conditions.

You also benefit from our **Erlab Safety Program** which includes **eValiQuest**, **ValiPass** and **ValiGuard**.

By subscribing to the Erlab offer, you gain access to the manufacturer's maintenance offer.

An offer customised for you

Professional installation, optimised filter replacement, compliance checks and scheduled maintenance are essential to ensure maximum user safety, long-lasting equipment performance and control over operating costs.

Our services (Only in France)

Installation and training

A certified technician will unpack and install your Erlab product in order to integrate it into your current installation and perform compliance checks. You will be trained on how to use it and its range of features, as well as in best practices so you can optimise the performance of your equipment.

Our contract offer (Only in France)

The annual maintenance includes an analysis of the evolution of your handlings to ensure that your configuration corresponds to your use and is in accordance with ValiPass (risk analysis certificate).

Inspection of your device is carried out according to a strict protocol, ensuring compliance with the NF X 15-211 standard.

This maintenance includes checking the air face velocity, the filtration efficiency, as well as the dynamic flow through a smoke test.



ERLAB FILTERS

Online shop for Erlab filters


Erlab has developed a **unique range of activated carbons** dedicated to inhalation protection for laboratory staff.


However, activated carbon is not the only guarantee of filtration performance. The efficiency of a molecular filter also depends on its installation and use.

The compliance of our hoods with the **AFNOR NF X 15-211** standard guarantees molecular filtration performance.

Purchase your filters online at filters.erlab.com

Please contact us for more information.

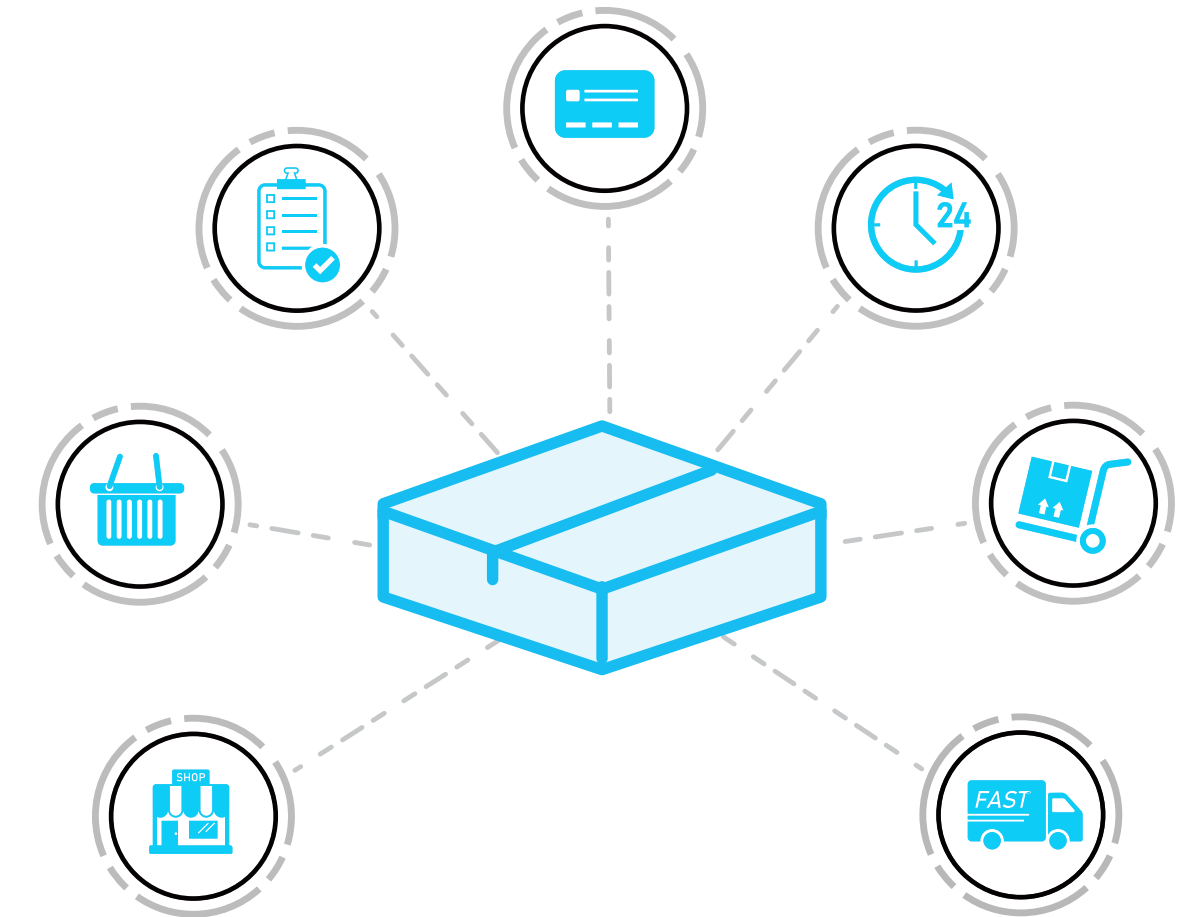

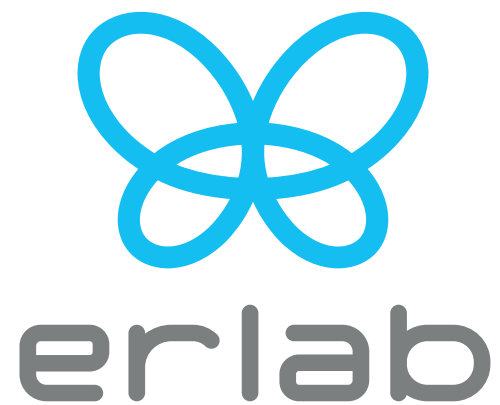

ERLAB D.F.S. S.A.S
Parc d'Affaires des Portes
B.P. 403 - 27104 Val-de-Reuil - France
www.erlab.com



France ventes@erlab.net

United States captairsales@erlab.com

China sales.china@erlab.com.cn

Germany export.north@erlab.net

United Kingdom export.north@erlab.net

Spain export.south@erlab.net

Italy export.south@erlab.net